



ペプチド - LC/MS
Peptides - LC/MS

A020808A

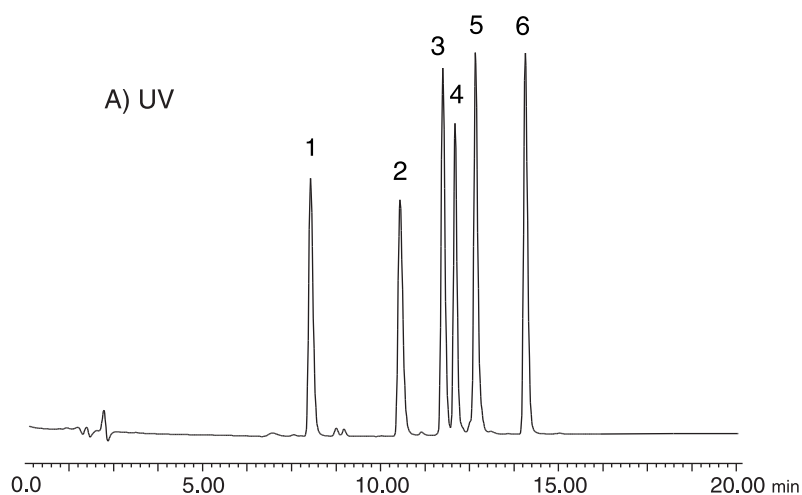

1. [Arg8]-Vasopressin
2. Bradykinin
3. Angiotensin II
4. Oxytocin
5. Angiotensin I
6. Leu-Enkephalin

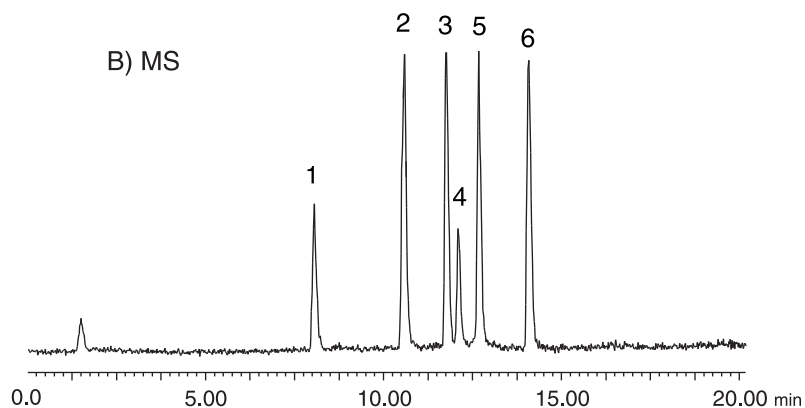

Column : YMC-Pack Pro C18 RS (5 μ m, 8nm)
150 \times 2.0 mm I.D.

Eluent : A) water/TFA (100/0.01)
B) acetonitrile/TFA (100/0.01)
10-35%B(0-15 min), 35%B(15-20 min)

Flow rate : 0.2 mL/min

Temperature : 37 $^{\circ}$ C

Detection : A) UV at 220 nm, B) ESI positive-mode

Injection : 4 μ L (0.05 mg/mL)