
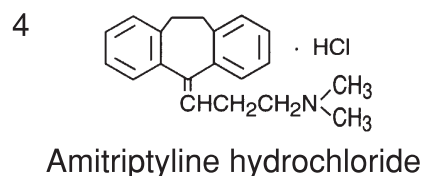
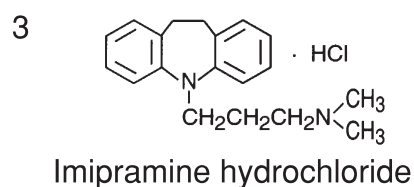
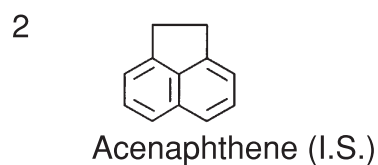
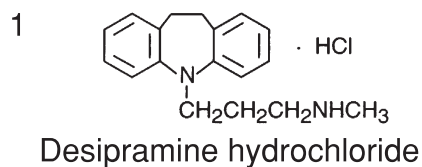


三環系抗うつ薬
Tricyclic antidepressants

A961016B

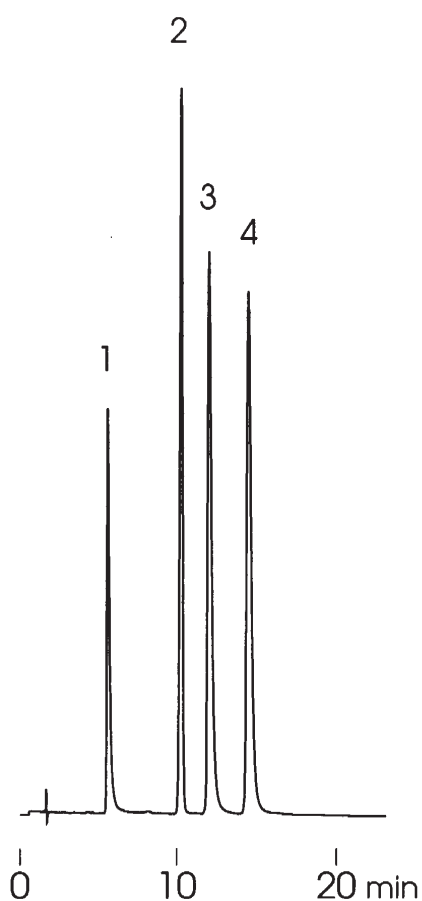

Column : YMC-Pack ODS-AM (5 μ m, 12nm)
150 \times 4.6 mm I.D.
Eluent : 20mM KH₂PO₄-K₂HPO₄ (pH 6.9)/methanol (25/75)
Flow rate : 1.0 mL/min
Temperature : 37 $^{\circ}$ C
Detection : UV at 254 nm, 0.16 AUFS
Injection : 5 μ L (0.1~0.5 mg/mL)