

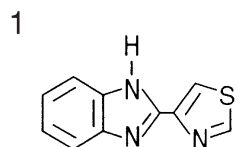


防かび剤  
Fungicides


A970217A

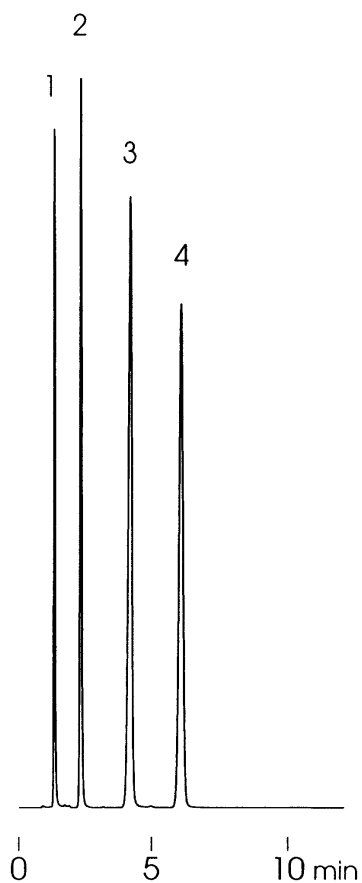


Thiabendazole

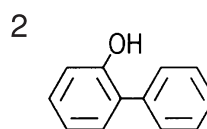

o-Phenylphenol

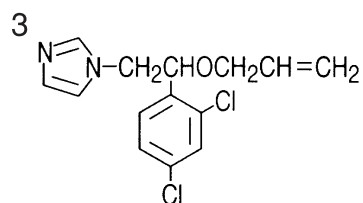

Imazalil

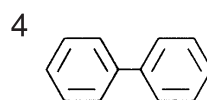

Diphenyl

Column : YMC-Pack Pro C18 (5  $\mu$ m, 12nm)  
75  $\times$  4.6 mm I.D.  
Eluent : methanol/20mM NH<sub>4</sub>H<sub>2</sub>PO<sub>4</sub> (70/30)  
Flow rate : 1.0 mL/min  
Temperature : 37  $^{\circ}$ C  
Detection : UV at 230 nm, 0.32 AUFS  
Injection : 10  $\mu$ L (0.02~0.20 mg/mL)