



利尿薬
Diuretics

A970426D

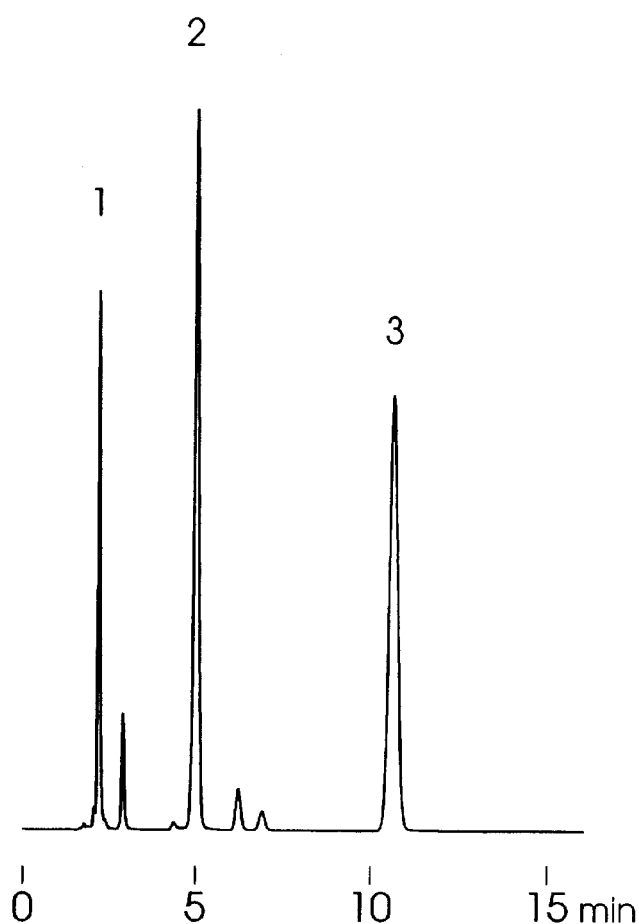


1

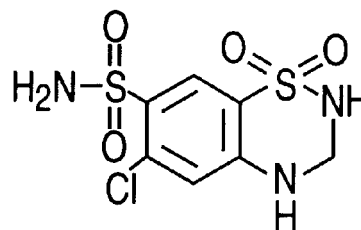


Hydrochlorothiazide

2

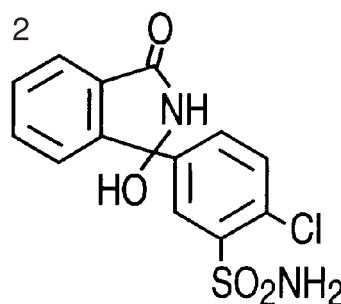

Chlorthalidone

3

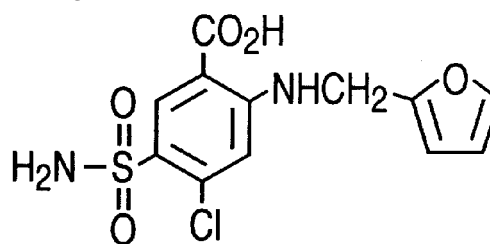

Furosemide

Column	: YMC-Pack ODS-AQ (5 μm, 12nm) 150 × 4.6 mm I.D.
Eluent	: 20mM NH ₄ H ₂ PO ₄ - H ₃ PO ₄ (pH 3.2)/methanol (60/40)
Flow rate	: 1.0 mL/min
Temperature	: 37 °C
Detection	: UV at 280 nm, 0.16 AUFS
Injection	: 15 μL (0.015~0.2 mg/mL)