


オリゴ糖
Oligosaccharides

A970701G

1. Fruf β 2 \rightarrow 1Fruf β 2 \rightarrow 1 α Glc (1-Kestose)
2. Fruf β 2 \rightarrow [1Fruf β 2]₂ \rightarrow 1 α Glc (Nystose)

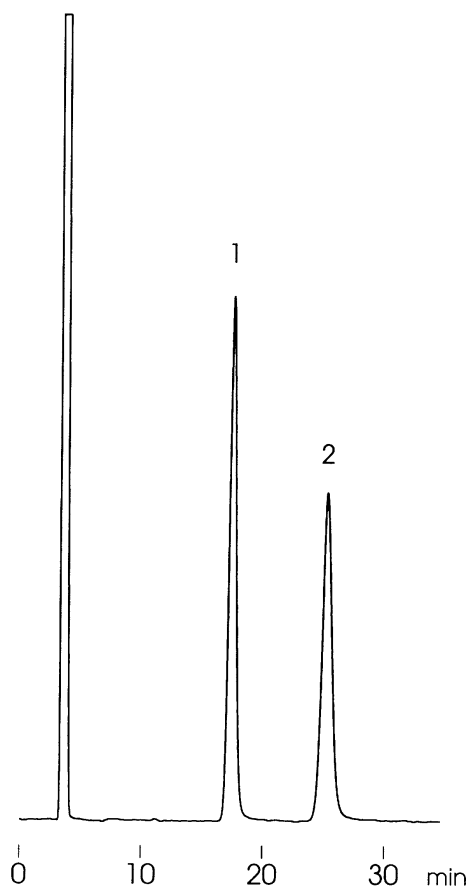

Column : YMC-Pack NH₂ (5 μ m, 12nm)
250 \times 4.6 mm I.D.
Eluent : acetonitrile/water (70/30)
Flow rate : 1.0 mL/min
Temperature : ambient (27 $^{\circ}$ C)
Detection : RI, 32 \times 10⁻⁶ RIU/FS
Injection : 20 μ L (4.0 mg/mL)