
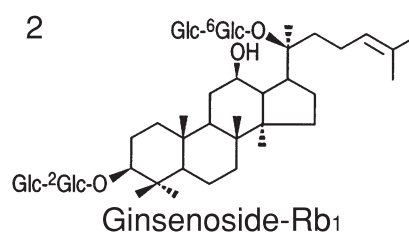
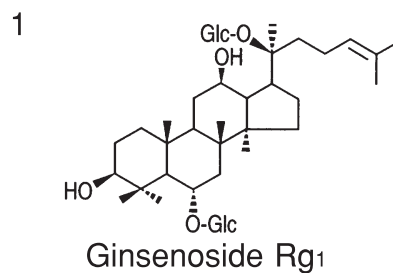


ジンセノサイド
Ginsenosides

A971121D

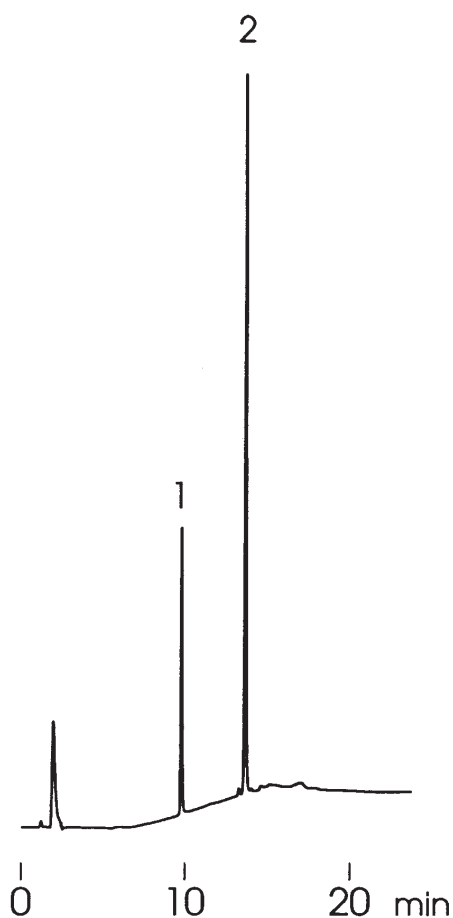

Column	: YMC-Pack Pro C18 (5 μ m, 12nm) 150 \times 3.0 mm I.D.
Eluent	: A)acetonitrile/water (20/80) B)acetonitrile/water (40/60) 0-100%B (0-10min), 100%B (10-20min)
Flow rate	: 0.4 mL/min
Temperature	: 37 $^{\circ}$ C
Detection	: UV at 203 nm, 0.32 AUFS
Injection	: 4 μ L (0.125, 0.5 mg/mL)