



バルサルタン・バルサルタン錠（日本薬局方収載原案記載条件）

Valsartan · Valsartan tablets (The draft for the Japanese Pharmacopoeia)

D120111A

Standard solution\*

(0.05 mg/mL Valsartan, 0.06 mg/mL Diclofenac sodium)

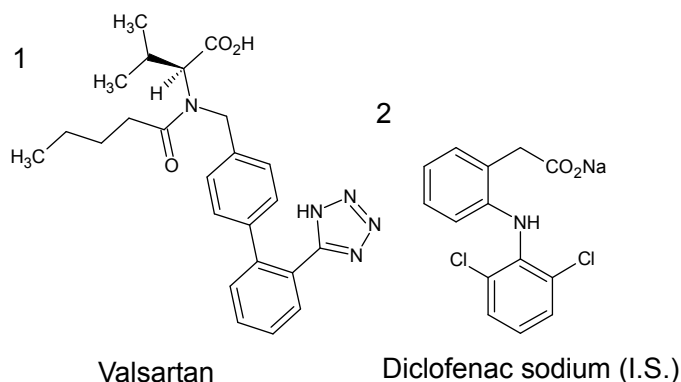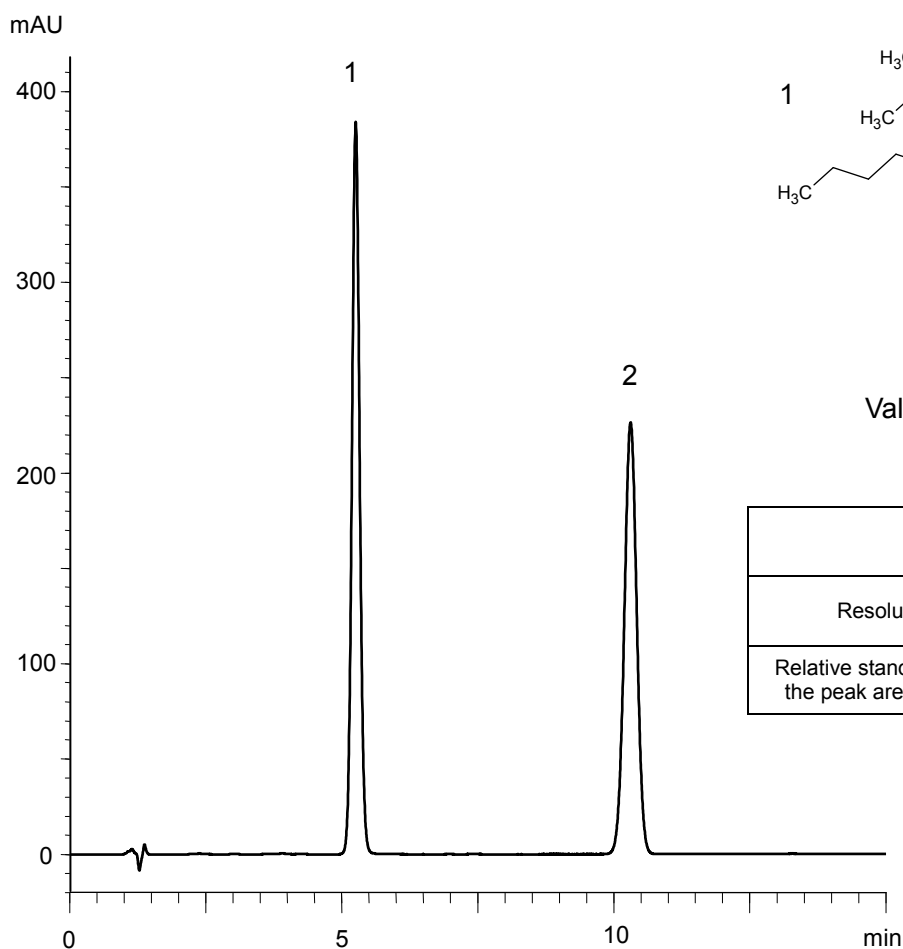

	System suitability requirement	Result
Resolution (1, 2)	$\geq 5$	14.9
Relative standard deviation of the peak area ratio of 1 to 2	$\leq 1.0\%$	0.05%

Column : YMC-Triart C18 (5  $\mu$ m, 12 nm)  
125 X 3.0 mm I.D.

Eluent : acetonitrile/water/acetic acid (50/50/0.1)

Flow rate : 0.45 mL/min (adjust the flow rate so that the retention time of valsartan is about 5 min)

Temperature : 25°C

Detection : UV at 225 nm

Injection : 10  $\mu$ L

(The draft for the Japanese Pharmacopoeia; Assay)

\* Standard solution was prepared from Valsartan supplied as a reagent for laboratory use.