
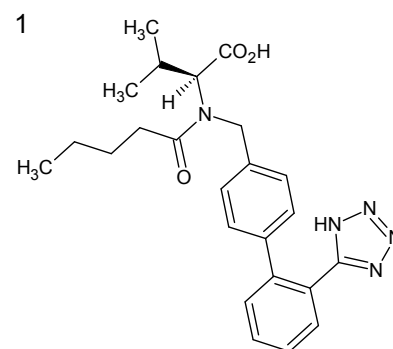



バルサルタン (米国薬局方記載条件)
Valsartan (The United States Pharmacopeia)

E130620B

Test solution*
(1.0 mg/mL Valsartan)

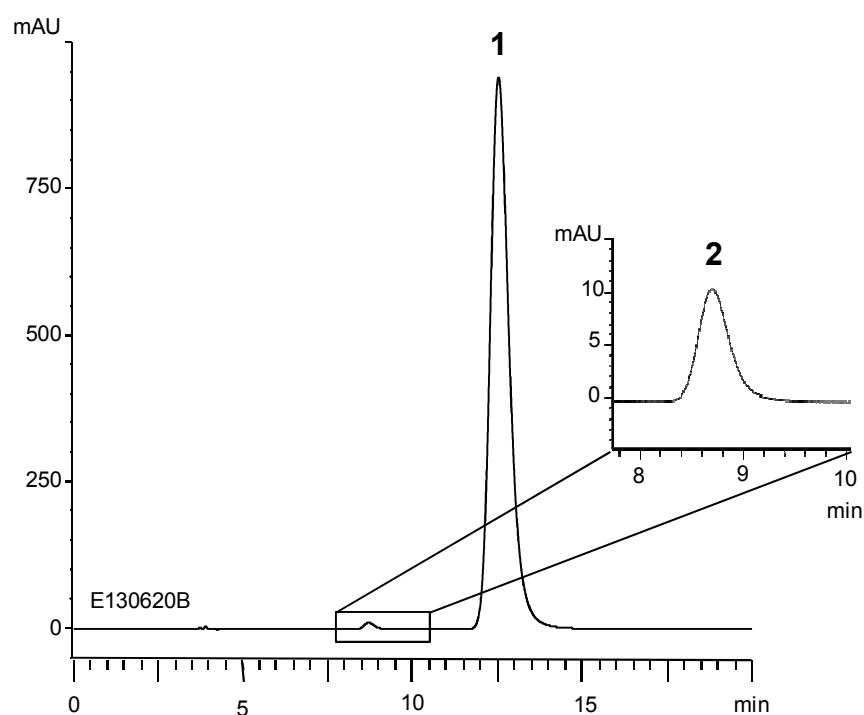

Valsartan

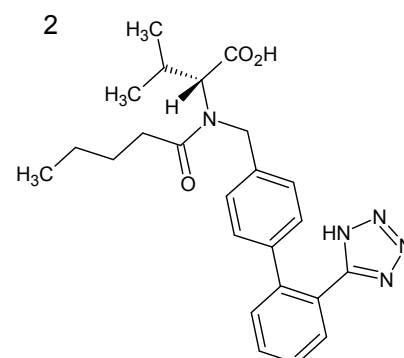

Valsartan enantiomer
(Valsartan related compound A)

Column : CHIRAL ART Cellulose-C (5 μ m)
250 X 4.6 mm I.D.
Eluent : *n*-hexane/2-propanol/trifluoroacetic acid (85/15/0.1)
Flow rate : 0.8 mL/min
Temperature : 25°C
Detection : UV at 230 nm
Injection : 10 μ L
(The United States Pharmacopeia 34th; Related compounds)

* Test solution was prepared from Valsartan supplied as a reagent for laboratory use.