

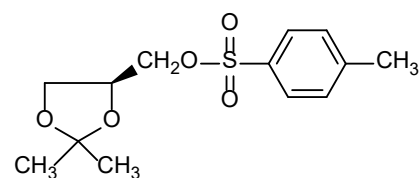


2,2-ジメチル-1,3-ジオキサラン-4-イルメチル *p*-トルエンスルホナート
2,2-Dimethyl-1,3-dioxolane-4-ylmethyl *p*-Toluenesulfonate

E180509F

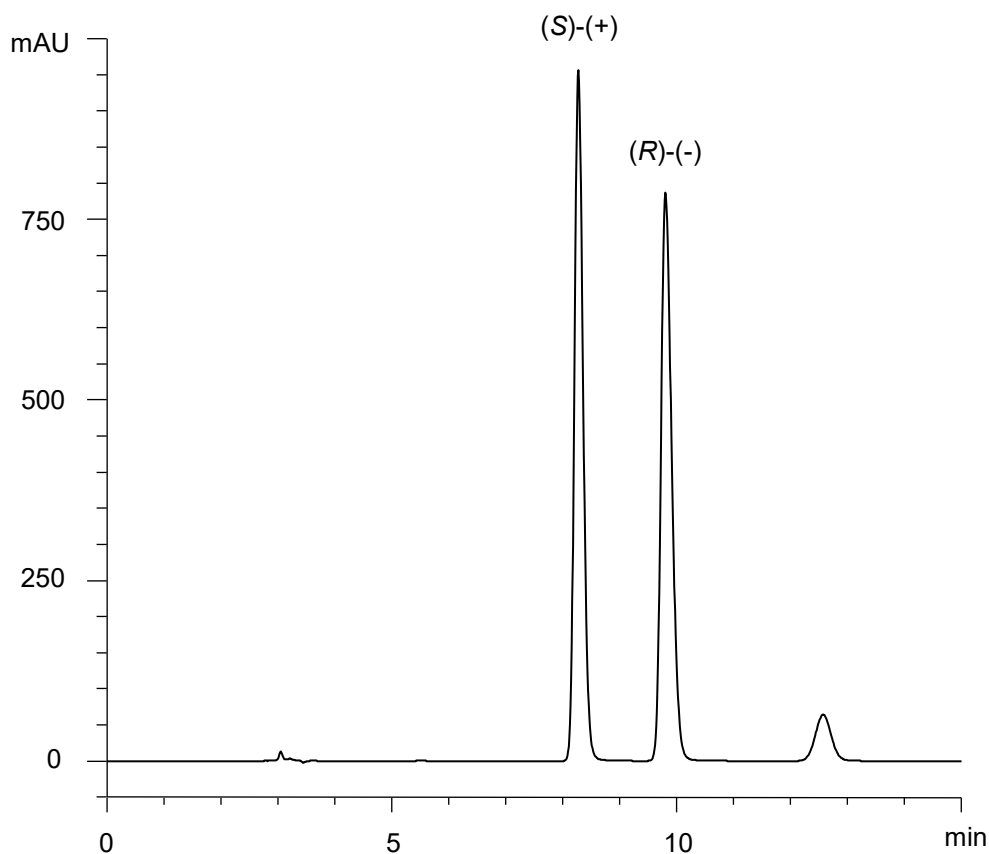

(S)-(+)

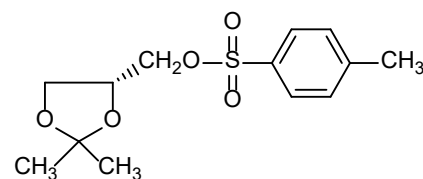

(R)-(-)

α	1.2
R_s	4.7

Column : CHIRAL ART Amylose-SA (5 μ m)
250 X 4.6 mm I.D.
Eluent : *n*-hexane/ethanol (80/20)
Flow rate : 1.0 mL/min
Temperature : 25°C
Detection : UV at 220 nm
Injection : 10 μ L (0.5 mg/mL)