


アルカロイド
Alkaloids

F090622C

1. Scopolamine
2. Atropine
3. Cinchonine
4. Quinine
5. Dihydroquinine

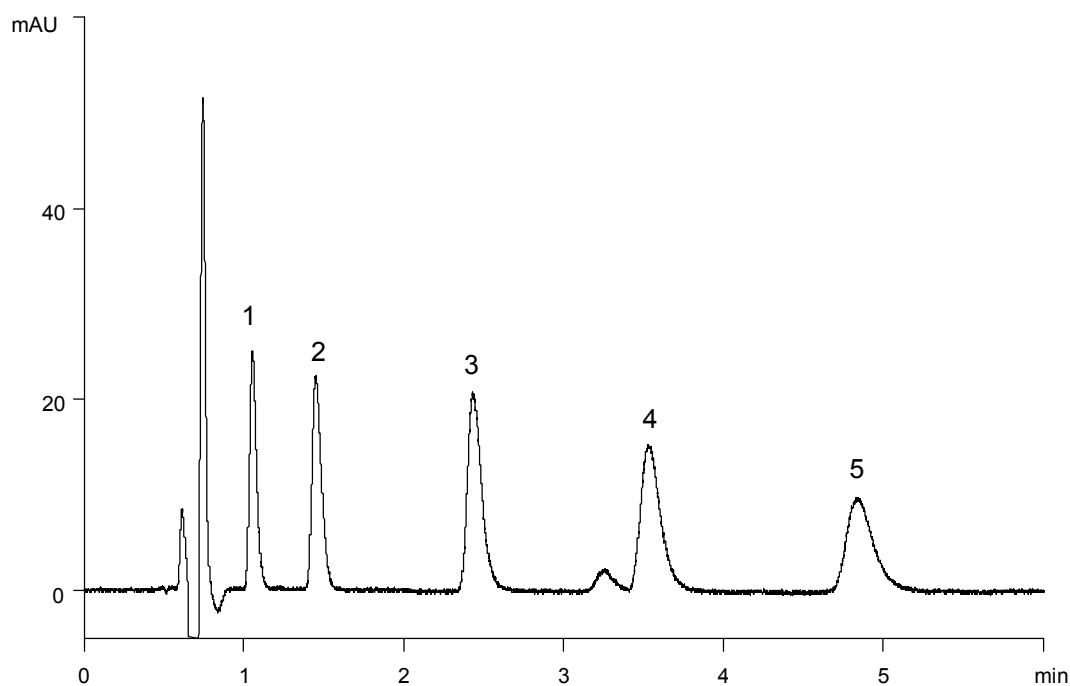

Column : YMC-Triart C18 (5 μ m, 12 nm)
50 X 2.0 mmI.D.

Eluent : 20 mM CH₃COOH-CH₃COONH₄ (pH 4.9)/acetonitrile (80/20)

Flow rate : 0.2 mL/min

Temperature : 40°C

Detection : UV at 220 nm

Injection : 1 μ L (0.02-0.1 mg/mL)