



ドネペジル塩酸塩錠・細粒（日本薬局方記載条件）
Donepezil hydrochloride tablets and fine granules
(The Japanese Pharmacopoeia)

F140228A

Standard solution*
(3.3 µg/mL Donepezil hydrochloride)

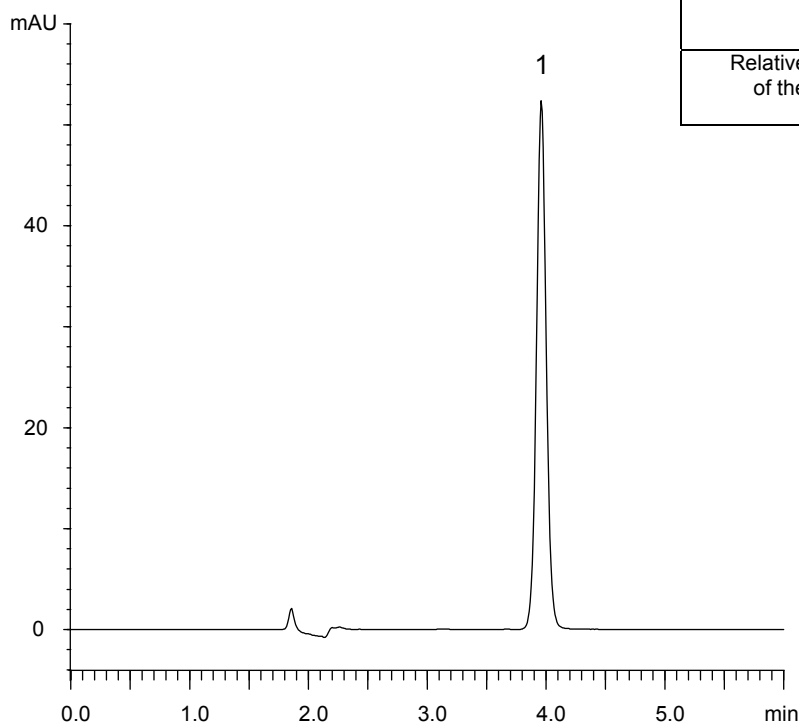

	System suitability requirement	Result
Theoretical plate number (Donepezil)	≥ 5000	9400
Tailing factor (Donepezil)	≤ 1.5	1.09
Relative standard deviation of the peak area (n=6) (Donepezil)	≤ 1.0 %	0.16 %

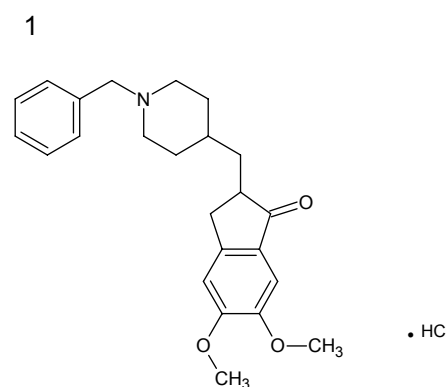

Donepezil hydrochloride

Column : YMC-Pack ODS-A (5 µm, 30 nm)
150 X 4.6 mm I.D.
Eluent : acetonitrile/water/perchloric acid (350/650/1)
Flow rate : 1.0 mL/min
(adjust the flow rate so that the retention time of donepezil is about 4 min)
Temperature : 35°C
Detection : UV at 271 nm
Injection : 50 µL
(The Japanese Pharmacopoeia 16th; Dissolution)

* Standard solution was prepared from Donepezil hydrochloride supplied as a reagent for laboratory use.