


インドメタシン
Indometacin

G900727F

1. *p*-Chlorobenzoic acid
2. *p*-Hydroxybenzoic acid n-butyl ester
3. Indometacin

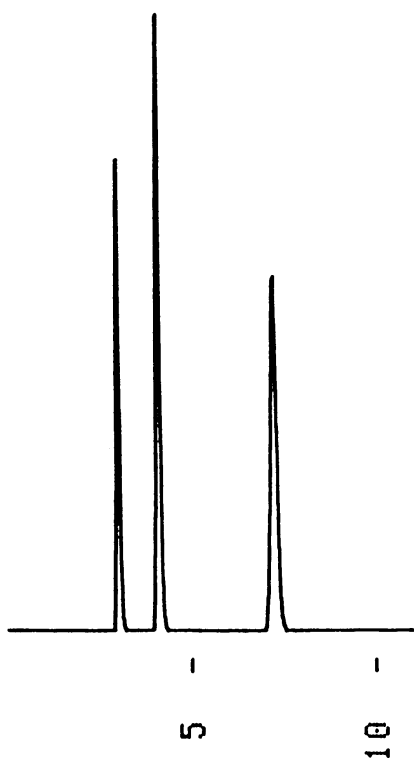

Column : YMC-Pack ODS-A (5 μ m, 12nm)
150 \times 4.6 mm I.D.
Eluent : methanol/0.1% phosphoric acid (70/30)
Flow rate : 1.3 mL/min
Temperature : 30 $^{\circ}$ C
Detection : UV at 254 nm, 0.64 AUFS
Injection : 20 μ L
(日本薬局方第12改正記載条件)