


インスリン
Insulins

G910522D

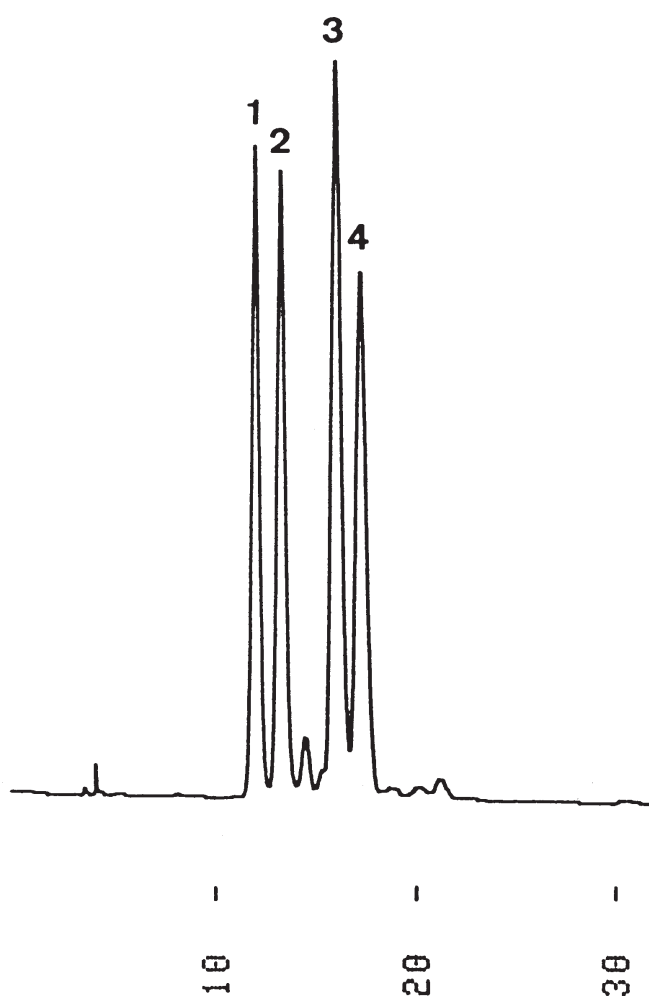

- 1. Insulin (Bovine)
- 2. Insulin (Sheep)
- 3. Insulin (Human)
- 4. Insulin (Porcine)

Column : YMC-Pack ODS-AM (5 μ m, 12nm)
150 \times 6.0 mm I.D.

Eluent : A) acetonitrile/0.01N HCl (28/72)
B) acetonitrile/0.01N HCl (35/65)
0-100%B (0-60min, linear)

Flow rate : 1.0 mL/min

Temperature : 30 $^{\circ}$ C

Detection : UV at 220 nm, 0.16 AUFS

Injection : 15 μ L (0.25 mg/mL)