


ペプチド
Peptides

G911217B

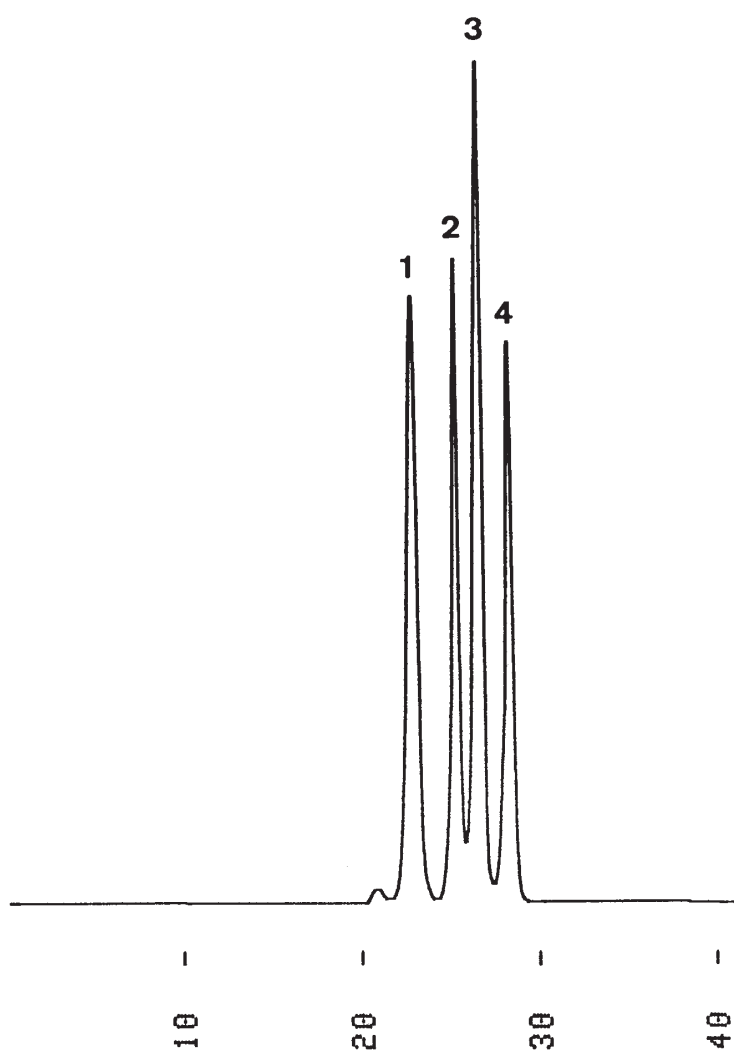

1. Insulin (Bovine)
2. Neurotensin
3. Angiotensin II
4. Glycine

Column	: YMC-Pack Diol-120 500 × 8.0 mm I.D.
Eluent	: 0.1M KH ₂ PO ₄ -K ₂ HPO ₄ (pH 7.0) containing 0.2M NaCl/acetonitrile (70/30)
Flow rate	: 0.7 mL/min
Temperature	: ambient (25 °C)
Detection	: UV at 215 nm, 0.16 AUFS
Injection	: 25 μL (0.07~5.3 mg/mL)