
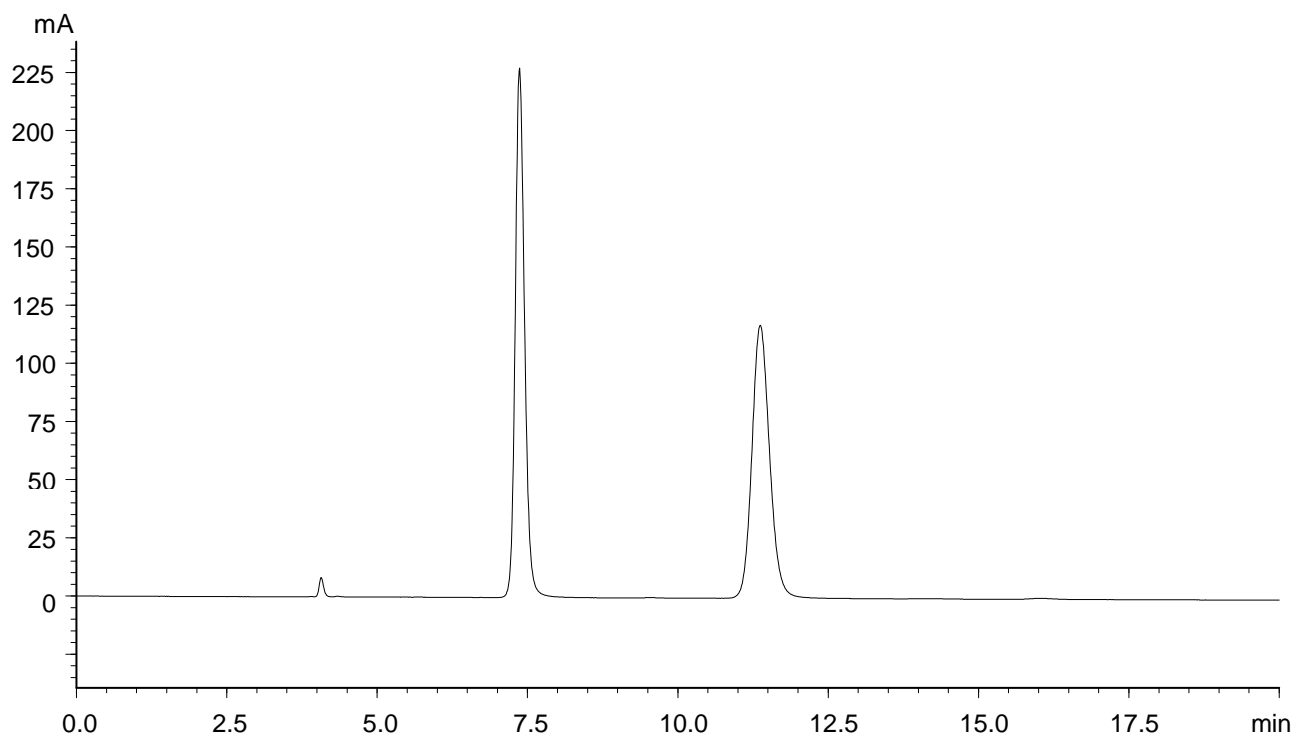


1,2-ビス[(2-メトキシフェニル)フェニルホスフィノ]エタン

1,2-Bis[(2-methoxyphenyl)phenyl-phosphino]ethane (DIPAMP)

H170821E

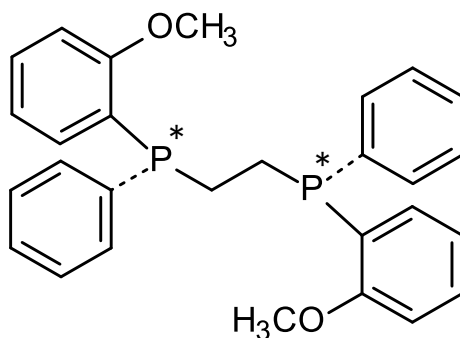

Column : CHIRAL ART Cellulose-SJ (5  $\mu$ m)  
250 X 4.6 mm I.D.

Eluent : *n*-hexane/chloroform (80/20)

Flow rate : 1.0 mL/min

Temperature : 25°C

Detection : UV at 290 nm

Injection : 5  $\mu$ L (1.0 mg/mL)