

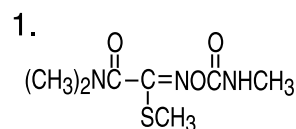


カルバメート系殺虫剤  
Carbamate insecticides


J010220D

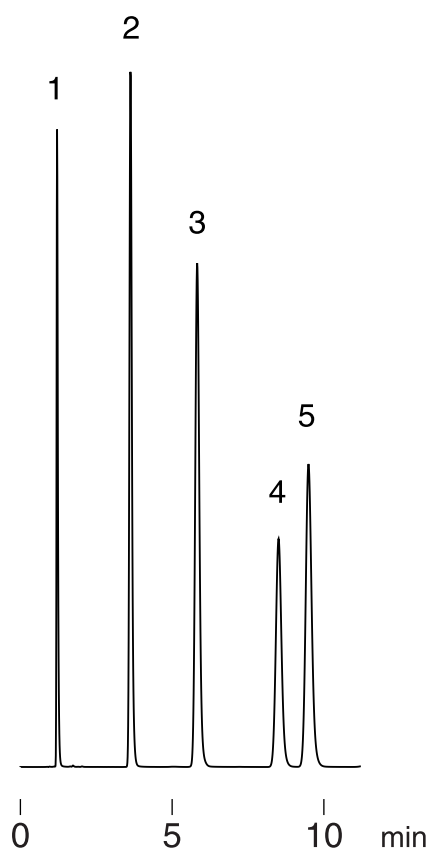


Oxamyl

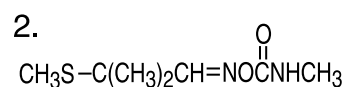


Aldicarb

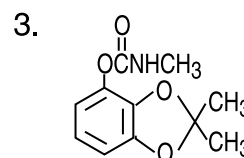

Bendiocarb

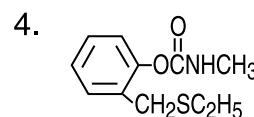

Ethiofencarb

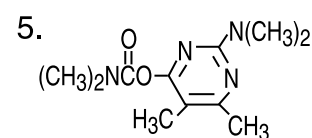

Pirimicarb

Column : Hydrosphere C18 (3  $\mu$ m, 12nm)  
75  $\times$  4.6 mm I.D.  
Eluent : methanol / 20mM KH<sub>2</sub>PO<sub>4</sub> (45/55)  
Flow rate : 1.0 mL/min  
Temperature : 30  $^{\circ}$ C  
Detection : UV at 240 nm, 0.64 AUFS  
Injection : 5  $\mu$ L (0.1~1.5 mg/mL)