
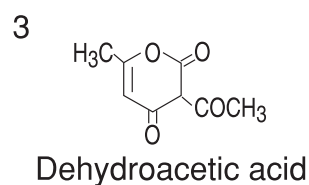
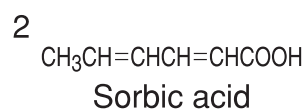
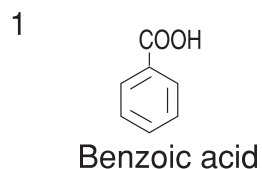


食品保存料  
Food preservatives

J011213F

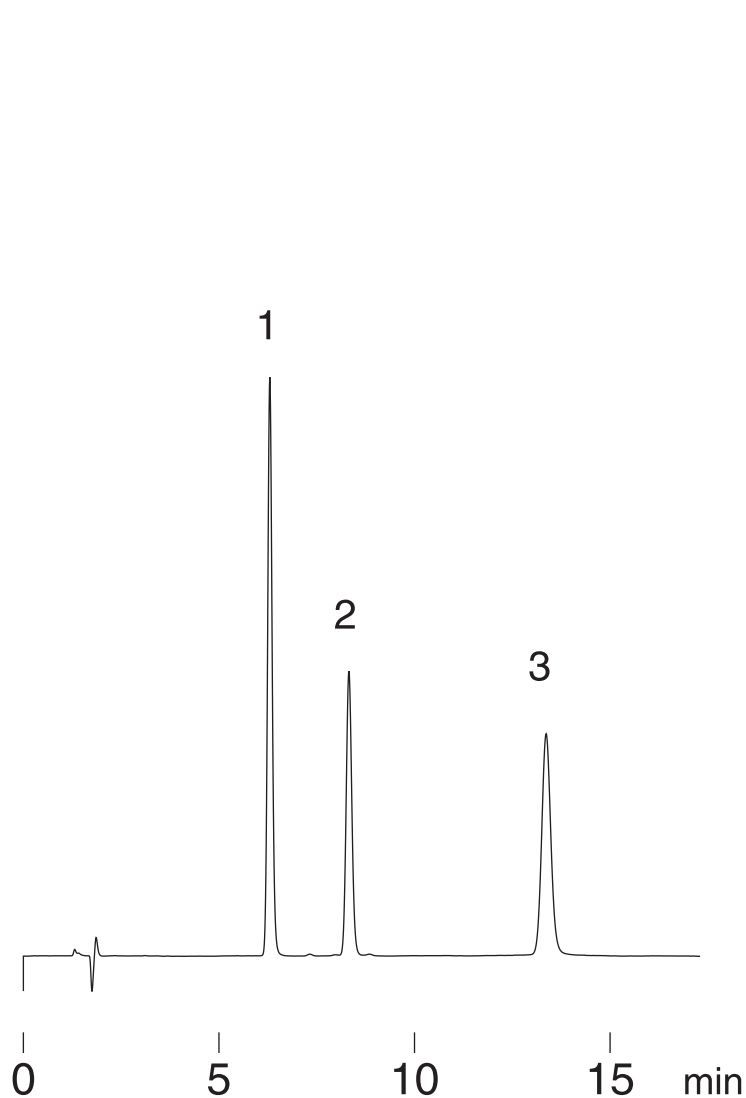

Column : YMC-Pack Pro C18 RS (5  $\mu$ m, 8nm)  
150  $\times$  4.6 mm I.D.  
Eluent : 20mM CH<sub>3</sub>COONa-CH<sub>3</sub>COOH (pH 4.3) / acetonitrile (80/20)  
Flow rate : 1.0 mL/min  
Temperature : 37  $^{\circ}$ C  
Detection : UV at 230 nm, 0.26 AUFS  
Injection : 10  $\mu$ L (0.02 mg/mL)