
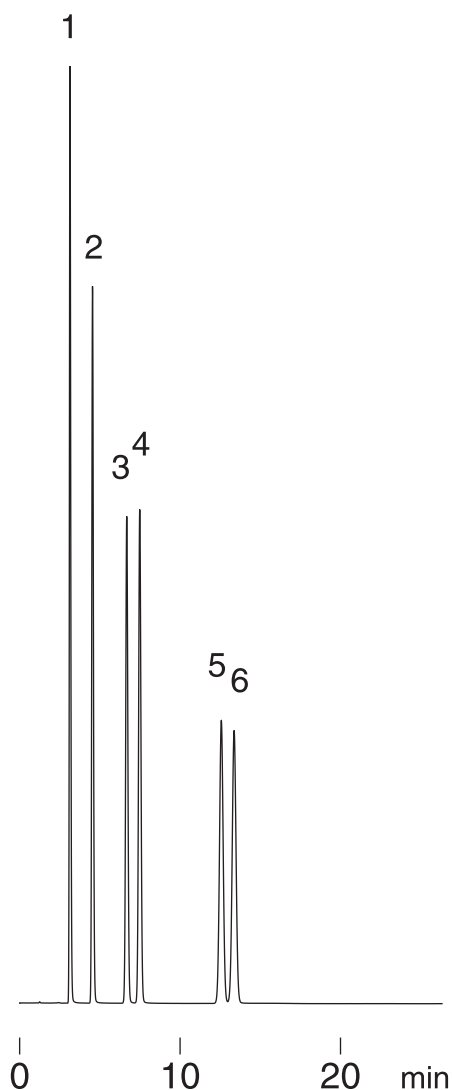
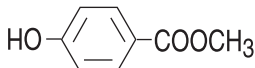
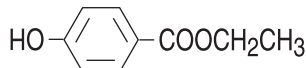
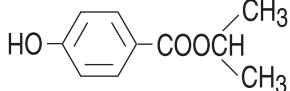
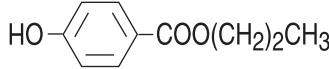
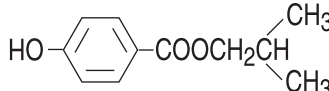


パラオキシ安息香酸エステル類
p-Hydroxybenzoic acid esters

J020826A

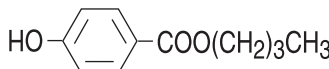

1. 
Methyl *p*-hydroxybenzoate
2. 
Ethyl *p*-hydroxybenzoate
3. 
Isopropyl *p*-hydroxybenzoate
4. 
Propyl *p*-hydroxybenzoate
5. 
Isobutyl *p*-hydroxybenzoate
6. 
Butyl *p*-hydroxybenzoate

Column : YMC-Pack Pro C18 RS (5 μ m, 8nm)
150 \times 4.6 mm I.D.
Eluent : methanol/water (60/40)
Flow rate : 1.0 mL/min
Temperature : 30 $^{\circ}$ C
Detection : UV at 270 nm, 0.32 AUFS
Injection : 4 μ L (0.1 mg/mL)