


合成抗菌剂
Antibacterial agents

K940614B

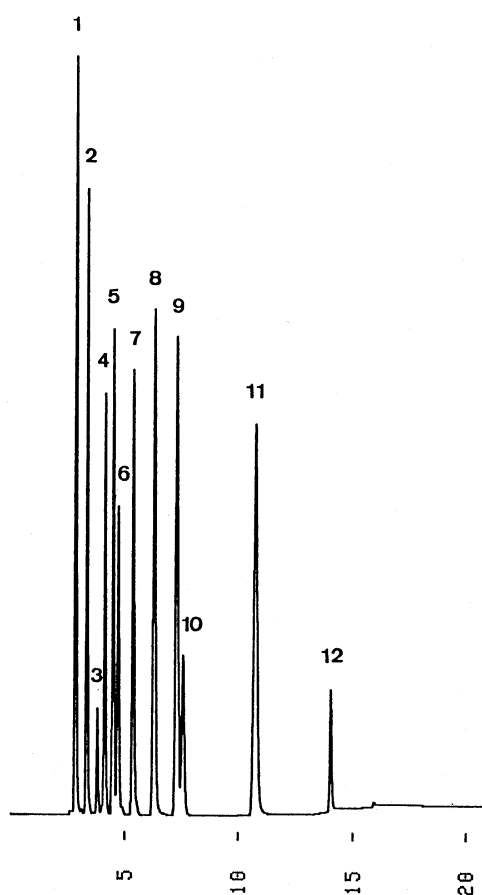

1. Olaquinox
2. Carbadox
3. Thiamphenicol
4. Sulfamerazine
5. Sulfadimidine
6. Sulfamonomethoxine
7. Furazolidone
8. Oxolinic acid
9. Sulfaquinoxaline
10. Sulfadimethoxine
11. Nalidixic acid
12. Piromidic acid

Column : YMC-Pack ODS-A (5 μ m, 12nm)
250 \times 4.6 mm I.D.

Eluent : A)acetonitrile/50mM NH₄H₂PO₄ (10/90)
B)acetonitrile/50mM NH₄H₂PO₄ (80/20)
25%B (0-5min), 25-100%B (5-15min), 100%B (15-20min)

Flow rate : 1.0 mL/min

Temperature : 37 $^{\circ}$ C

Detection : UV at 240 nm, 0.32 AUFS

Injection : 10 μ L (0.05 mg/mL)