


抗体薬物複合体
Antibody-Drug Conjugates

M180803EA

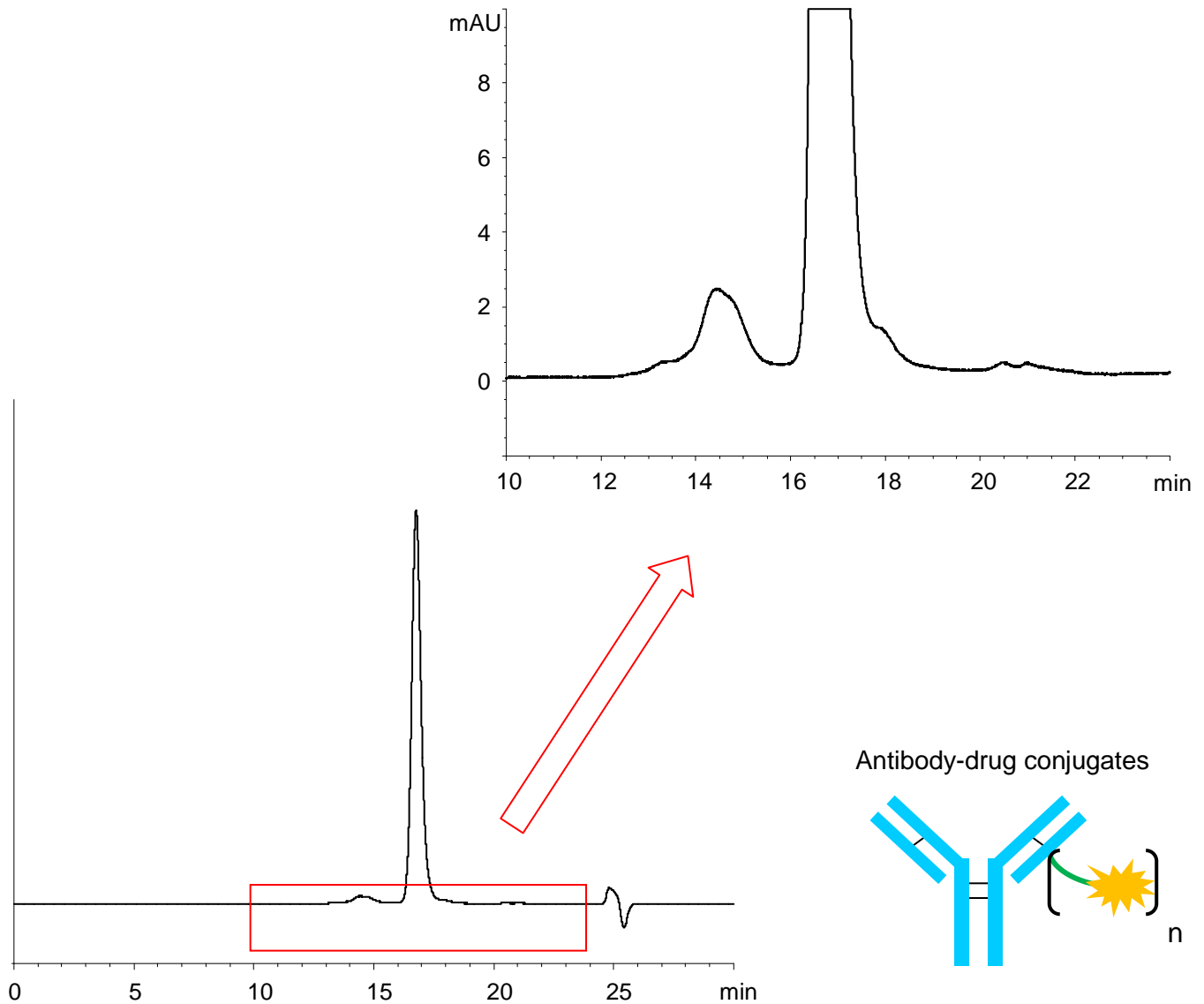

Column	: YMC-SEC MAB (3 μ m, 25 nm) 300 X 4.6 mm I.D.
Eluent	: 0.1 M KH_2PO_4 - K_2HPO_4 (pH 7.0) containing 0.2 M NaCl/2-propanol (85/15)
Flow rate	: 0.165 mL/min
Temperature	: 25°C
Detection	: UV at 280 nm
Injection	: 4 μ L (2.5 mg/mL)
Sample	: SigmaMAb Antibody Drug Conjugate (ADC) Mimic