


アミノアルコール
Aminoalcohols

M920124B

1. Ethanolamine
2. Diethanolamine
3. N-Methylethanolamine
4. Triethanolamine
5. N,N-Dimethylethanolamine

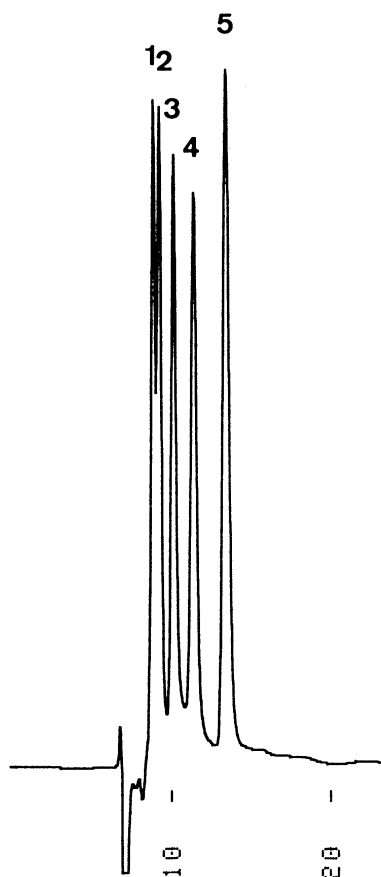

Column	: YMC-Pack PolymerC18 250 × 6.0 mm I.D.
Eluent	: 100mM Na ₂ HPO ₄ /100mM NaOH (60/40, pH 12.0)
Flow rate	: 0.6 mL/min
Temperature	: 20 °C
Detection	: UV at 215 nm, 0.32 AUFS
Injection	: 50 μL (0.2~3.0 mg/mL)