


水溶性ビタミン  
Water-soluble vitamins

M920207B

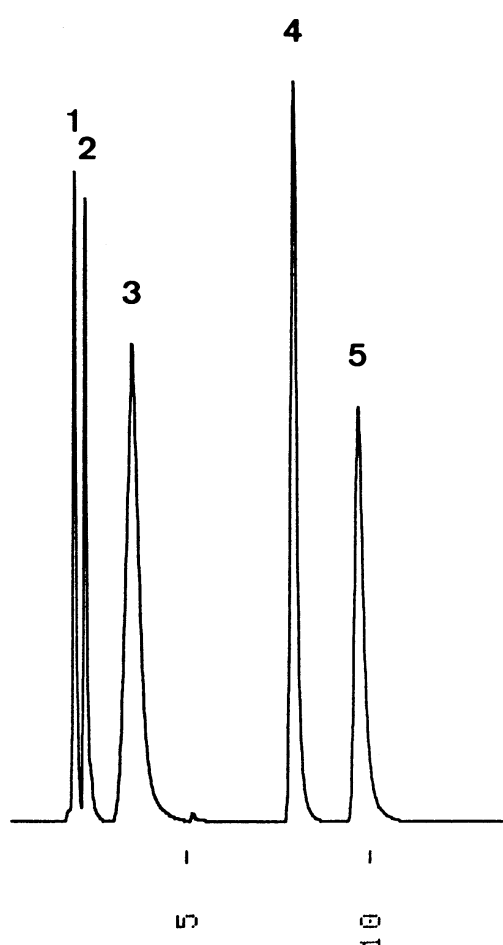

1. L-Ascorbic acid (Vitamin C)
2. Orotic acid (Vitamin B13)
3. Pyrroloquinoline quinone(PQQ)
4. Pyridoxine HCl (Vitamin B6)
5. Nicotinamide

Column : YMC-Pack ODS-AM (5  $\mu$ m, 12nm)  
150  $\times$  4.6 mm I.D.  
Eluent : 100mM CH<sub>3</sub>COOH/100mM CH<sub>3</sub>COONH<sub>4</sub> (30/70, pH 5.1)  
Flow rate : 1.0 mL/min  
Temperature : 37  $^{\circ}$ C  
Detection : UV at 254 nm, 0.16 AUFS  
Injection : 10  $\mu$ L (0.04~0.20 mg/mL)