


ペプチド(MW 556 - 3,465)

Peptides(MW 556 - 3,465)

N041021C

1. Oxytocin (MW 1,007)
2. Met-Enkephalin (MW 574)
3. Neurotensin (MW 1,673)
4. Leu-Enkephalin (MW 556)
5. γ -Endorphin (MW 1,859)
6. β -Endorphin(Human) (MW 3,465)

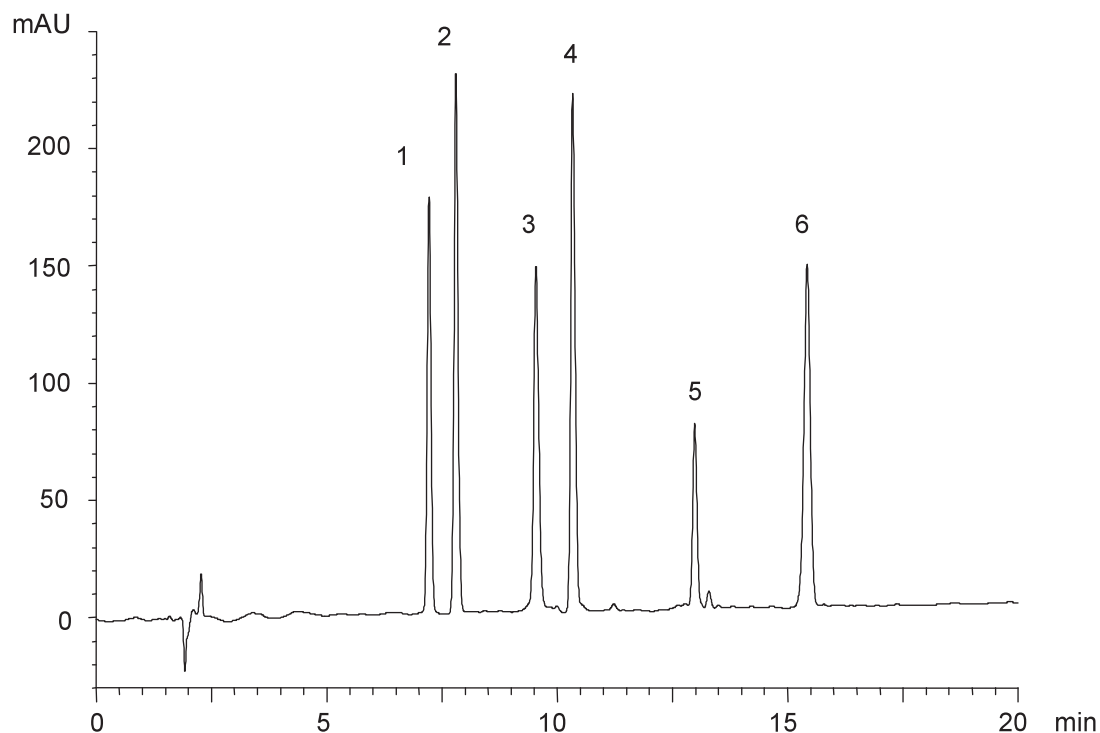

Column : YMC-Pack Pro C18 (5 μ m, 12nm)
150 \times 4.6 mm I.D.

Eluent : A) water / TFA (100/0.05)
B) acetonitrile / TFA (100/0.05)
20-40%B (0-20min)

Flow rate : 1.0 mL/min

Temperature : 37 $^{\circ}$ C

Detection : UV at 220 nm

Injection : 10 μ L (0.075~0.25 mg/mL)