
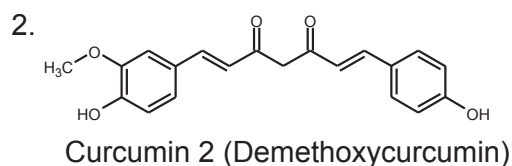
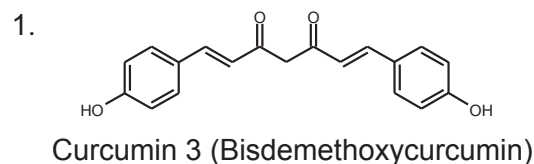


ウコンサプリメント中のクルクミン

Curcuminoids in a commercial turmeric supplement



N070729S

A) UV at 250 nm

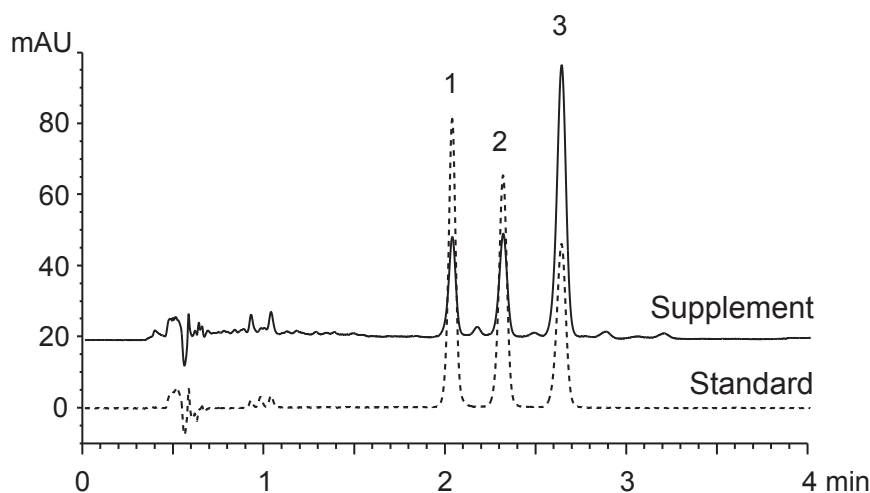

3

B) UV at 420 nm

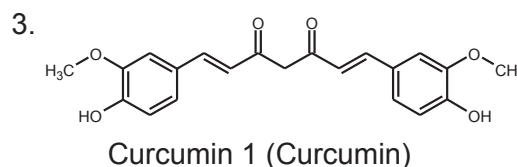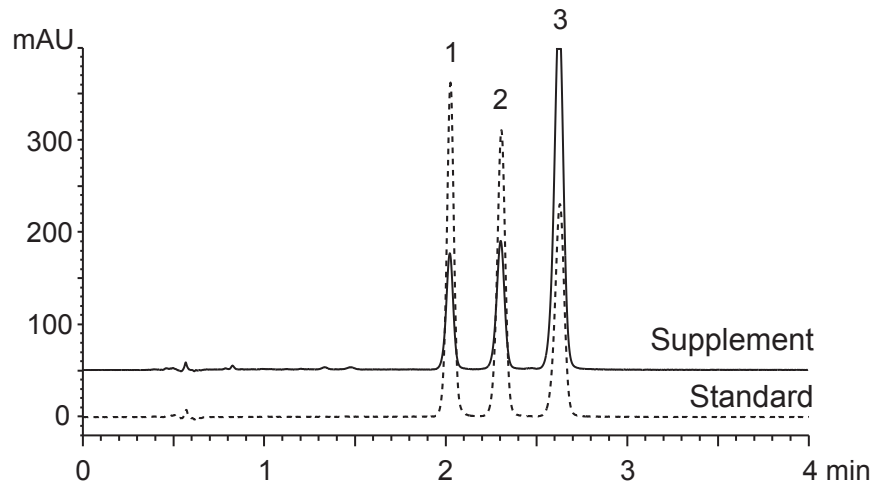

| | |
|-------------|---|
| Column | : YMC-UltraHT Pro C18 (2 μm, 12 nm) |
| Eluent | : 100 X 2.0 mm I.D. acetonitrile / water / formic acid (50/50/0.1) |
| Flow rate | : 0.4 mL/min |
| Temperature | : 40°C |
| Detection | : A) UV at 250 nm, B) UV at 420 nm |
| Injection | : standard 1 μL (0.1 mg/mL) : supplement 1 μL (6.9 mg/mL) |
| Sample | : methanol extract of a commercially available turmeric supplement |