
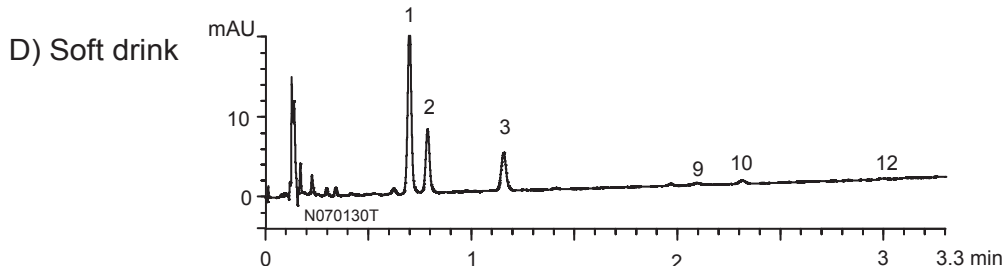
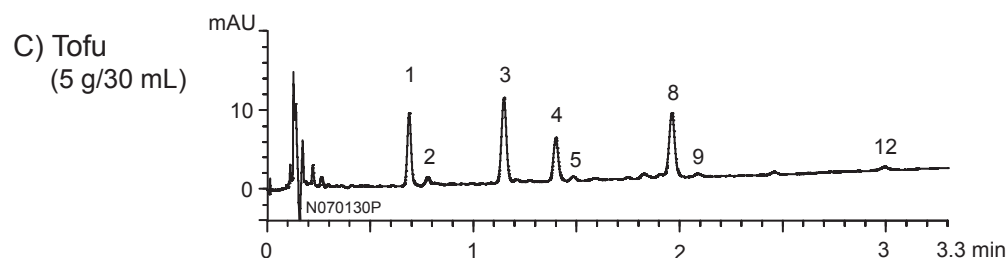
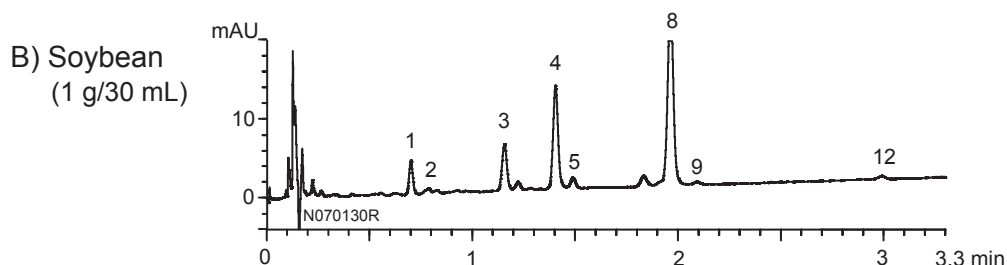



食品中の大豆イソフラボン  
Soy isoflavones in foods

N080415Z

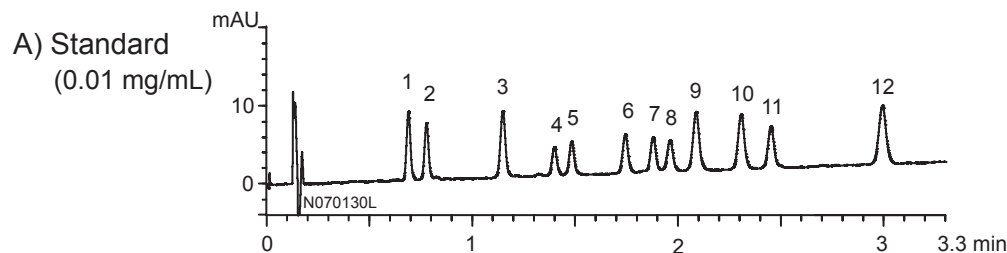

1. Daidzin
2. Glycitin
3. Genistin
4. 6''-O-Malonyldaidzin
5. 6''-O-Malonylglycitin
6. 6''-O-Acetyldaidzin
7. 6''-O-Acetylglycitin
8. 6''-O-Malonylgenistin
9. Daidzein
10. Glycitein
11. 6''-O-Acetylgenistin
12. Genistein

**Sample preparation method**

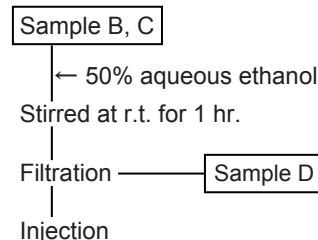

Column	: YMC-UltraHT Hydrosphere C18 (2 μm, 12 nm) 50 X 2.0 mmI.D.
Eluent	: A) water / acetic acid (100/3) B) acetonitrile / water / acetic acid (50/50/3) 25-60%B (0-3.3 min)
Flow rate	: 0.9 mL/min
Temperature	: 35°C
Detection	: UV at 254 nm
Injection	: 0.6 μL