


血漿成分
Plasma constituents

N080702A

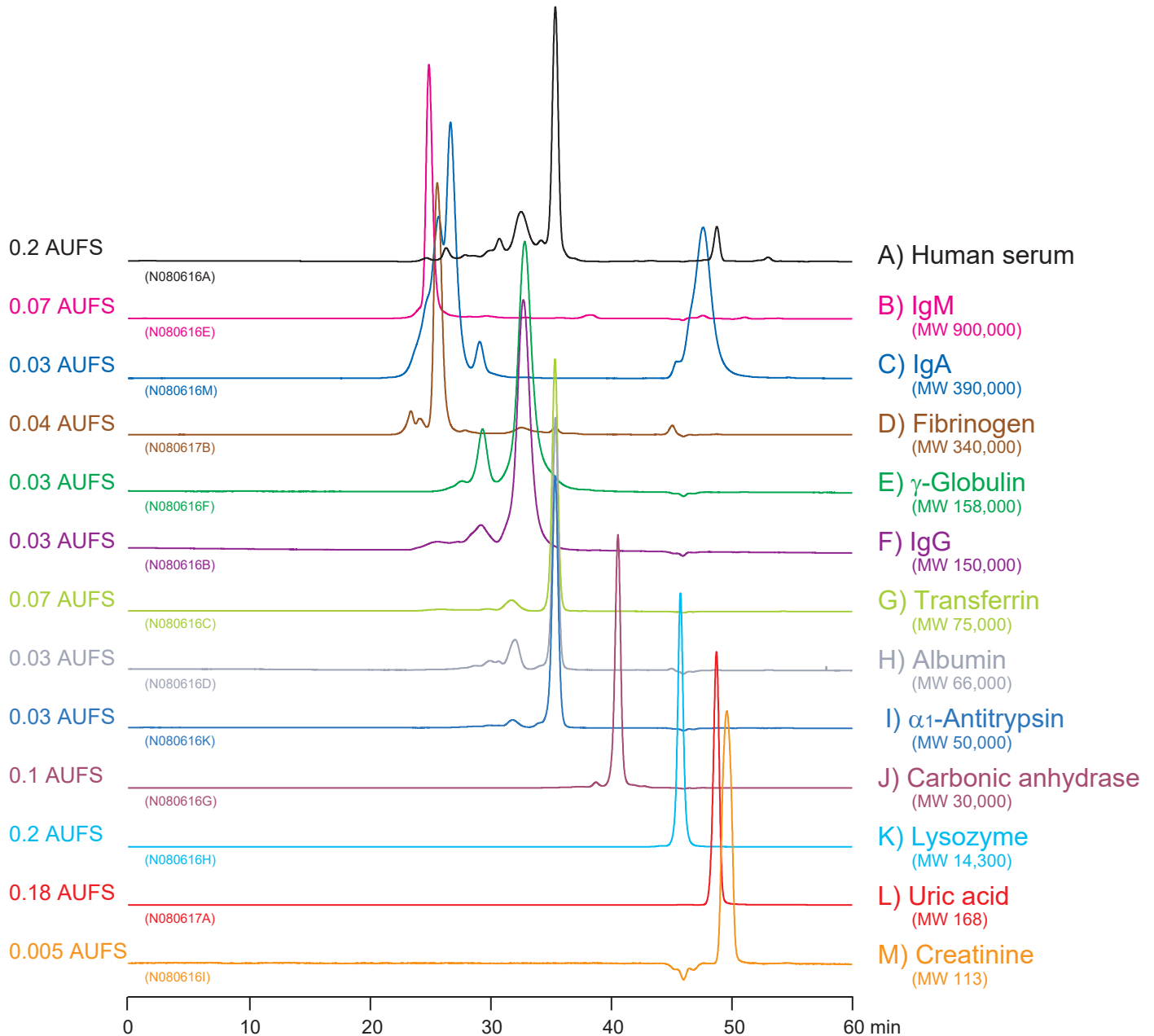

Column : YMC-Pack Diol-300 + Diol-200 (5 μ m)
 300 X 8.0 mm I.D. X 2
 Eluent : 0.1 M KH₂PO₄-K₂HPO₄ (pH 7.0) containing 0.2 M NaCl
 Flow rate : 0.5 mL/min
 Temperature : ambient (25°C)
 Detection : UV at 280 nm
 Injection : 20 μ L (L : 1 μ L)
 Sample : A) 100 μ L/mL
 B-M) 1.0 mg/mL