

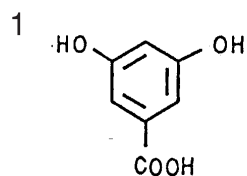


レゾルシン酸異性体
Resorcylic acid isomers


O940606J

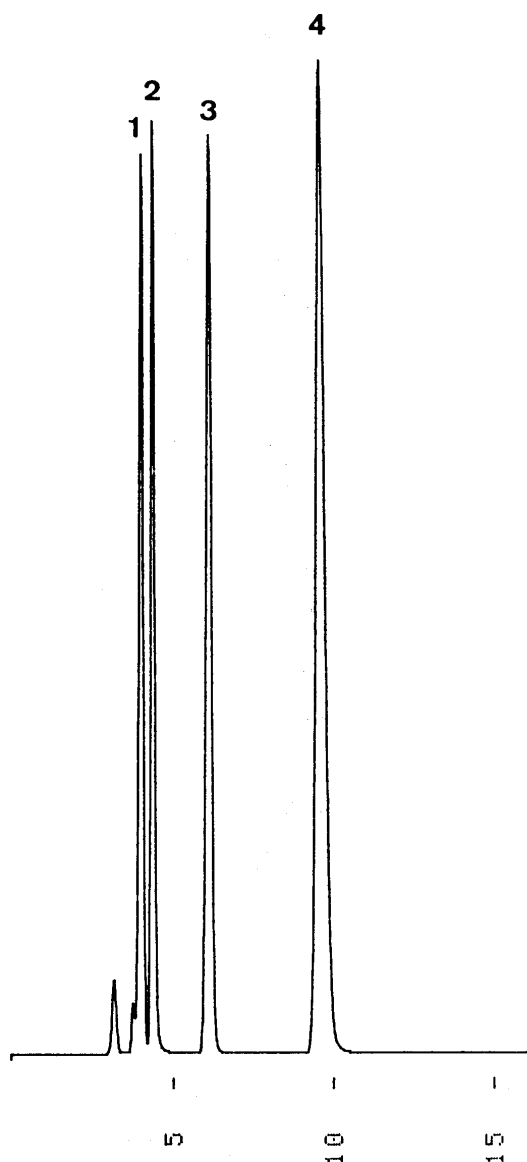


α -Resorcylic acid

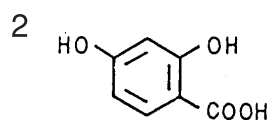

β -Resorcylic acid

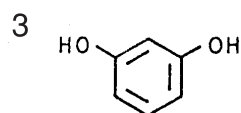

Resorcinol

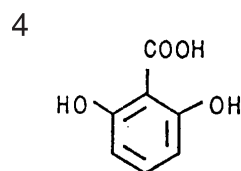

γ -Resorcylic acid

Column : YMC-Pack Ph (5 μ m, 12nm)
150 \times 4.6 mm I.D.
Eluent : methanol/100mM KH₂PO₄ (10/90)
Flow rate : 0.8 mL/min
Temperature : 37 $^{\circ}$ C
Detection : UV at 285 nm, 0.16 AUFS
Injection : 10 μ L (0.05~0.4 mg/mL)