


セフェム系抗生物質  
Cephalosporin antibiotics

O940718B

1. Cefsulodin sodium
2. Cephalexin hydrate
3. Cefaclor
4. Cephradine
5. Cephaloglycin
6. Cephapirin sodium
7. Cefotaxime sodium
8. Cephaloridine
9. Cefazolin sodium
10. Cefuroxime sodium
11. Cefoxitin sodium

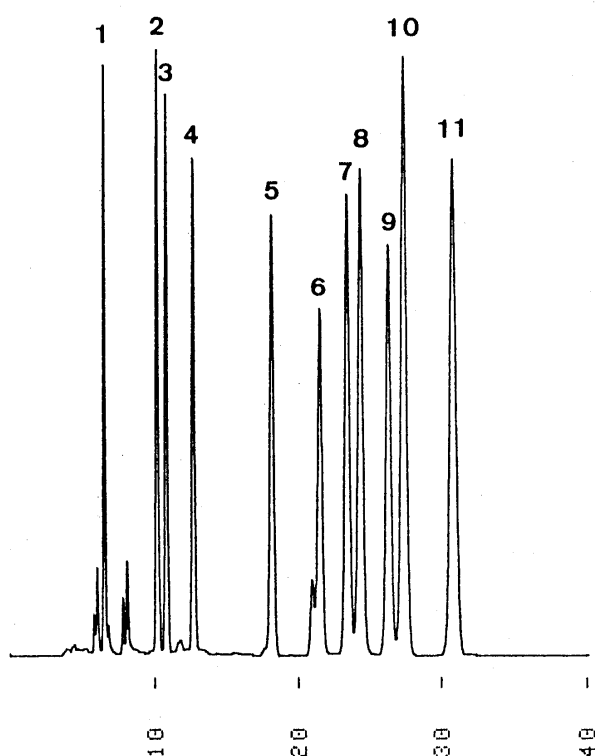

Column : J'sphere ODS-L80 (4  $\mu$ m, 8nm)  
250  $\times$  4.6 mm I.D.  
Eluent : acetonitrile/100mM KH<sub>2</sub>PO<sub>4</sub> (10/90)  
Flow rate : 0.8 mL/min  
Temperature : 37  $^{\circ}$ C  
Detection : UV at 260 nm, 0.16 AUFS  
Injection : 5  $\mu$ L (0.12~0.35 mg/mL)