

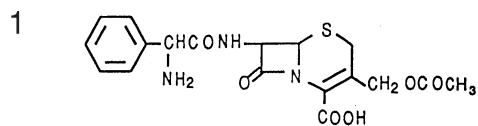


セフェム系抗生物質
Cephalosporin antibiotics


O940720A

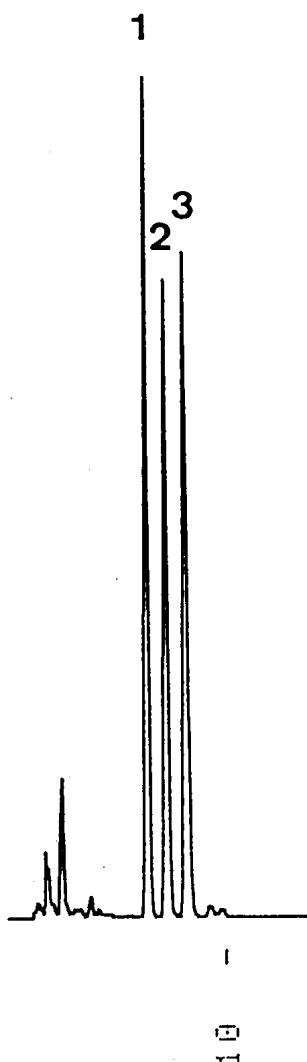

Cephaloglycin

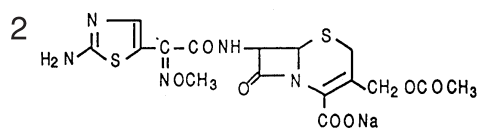

Cefotaxime sodium

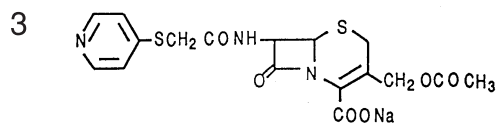

Cephapirin sodium

Column : J'sphere ODS-M80 (4 μ m, 8nm)
75 \times 4.6 mm I.D.
Eluent : acetonitrile/100mM KH₂PO₄ (10/90)
Flow rate : 0.8 mL/min
Temperature : 37 $^{\circ}$ C
Detection : UV at 260 nm, 0.64 AUFS
Injection : 5 μ L (0.5 mg/mL)