
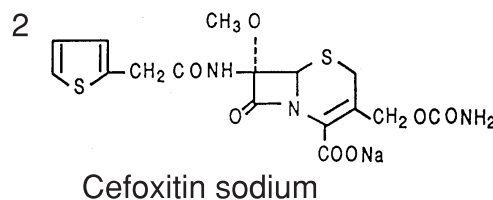
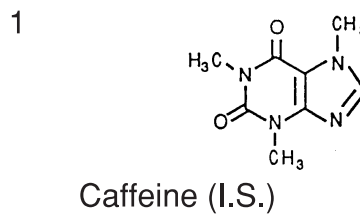


セフォキシチンナトリウム  
Cefoxitin sodium

O940729A

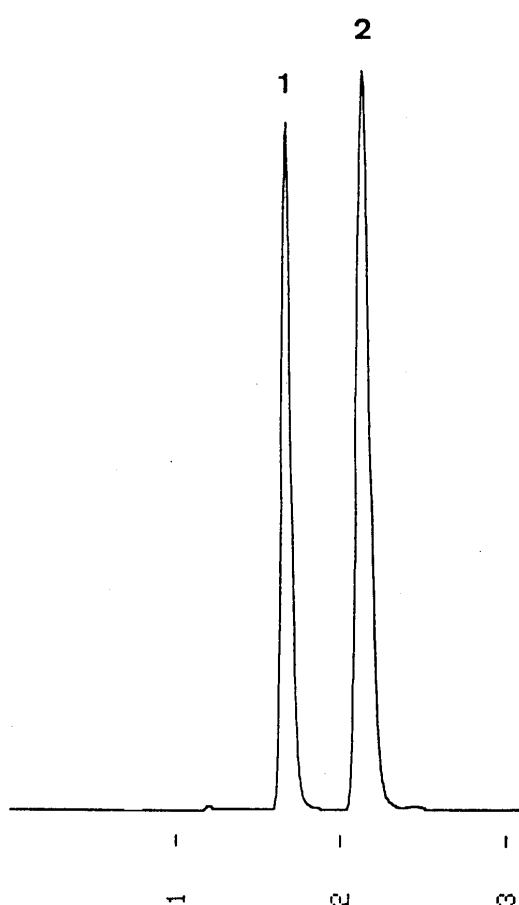

Column : J'sphere ODS-H80 (4  $\mu$ m, 8nm)  
150  $\times$  2.0 mm I.D.  
Eluent : acetonitrile/50mM KH<sub>2</sub>PO<sub>4</sub> (20/80)  
Flow rate : 0.34 mL/min  
Temperature : 37  $^{\circ}$ C  
Detection : UV at 260 nm, 0.32 AUFS  
Injection : 1.2  $\mu$ L (1.6, 2.4 mg/mL)