


ペプチド
Peptides

O940804B

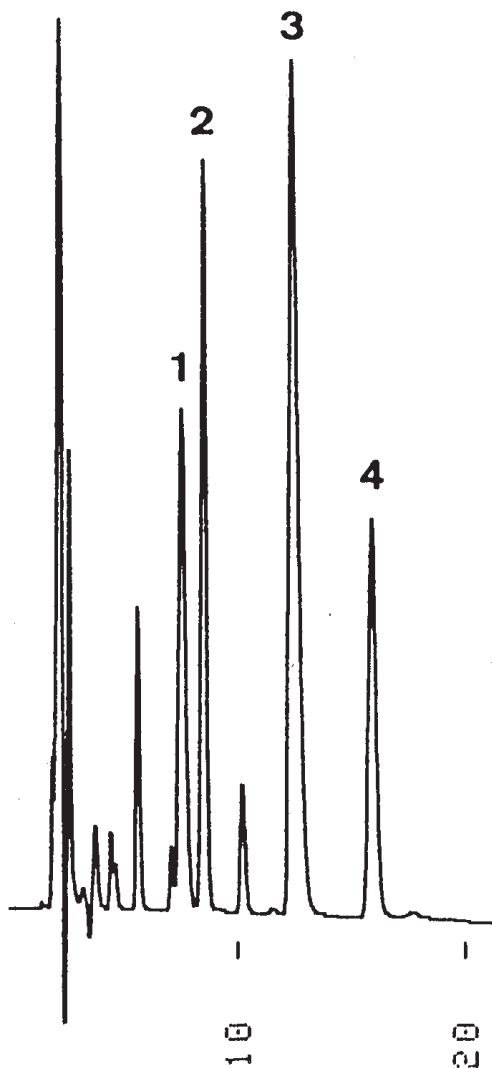

1. Bradykinin
2. Met-Enkephalin
3. Angiotensin III
4. Leu-Enkephalin

Column : J'sphere ODS-H80 (4 μ m, 8nm)
150 \times 1.5 mm I.D.

Eluent : acetonitrile/water/TFA (20/80/0.1)

Flow rate : 0.1 mL/min

Temperature : 37 $^{\circ}$ C

Detection : UV at 220 nm, 0.08 AUFS

Injection : 2 μ L (34~65 μ g/mL)