

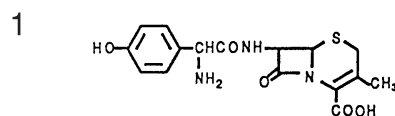


セフェム系抗生物質
Cephalosporin antibiotics


O940927B

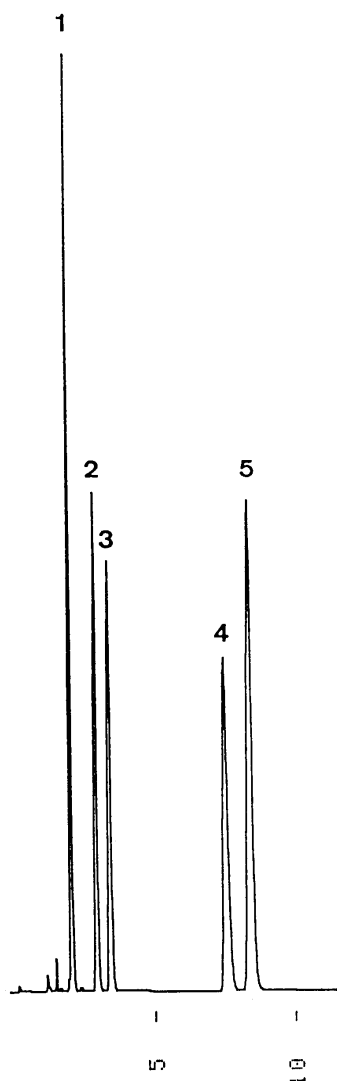


Cefadroxil

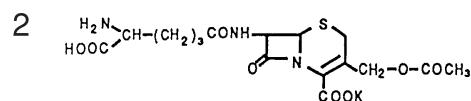


Cephalosporin C potassium

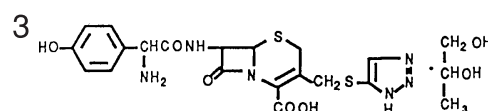

Cefatrizine propylene glycol

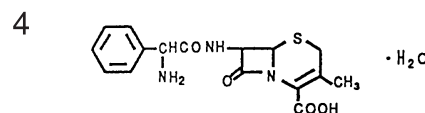

Cephalexin hydrate

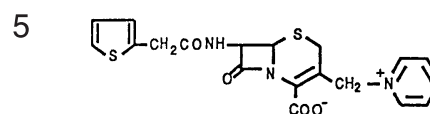

Cephaloridine

Column : J'sphere ODS-H80 (4 μ m, 8nm)
150 \times 4.6 mm I.D.
Eluent : methanol/water/acetic acid (10/85/5)
Flow rate : 1.0 mL/min
Temperature : 37 $^{\circ}$ C
Detection : UV at 260 nm, 0.16 AUFS
Injection : 6 μ L (0.2 mg/mL)