

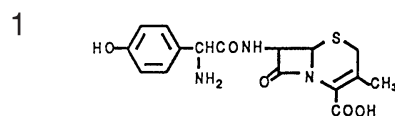


セフェム系抗生物質  
Cephalosporin antibiotics


O940927I

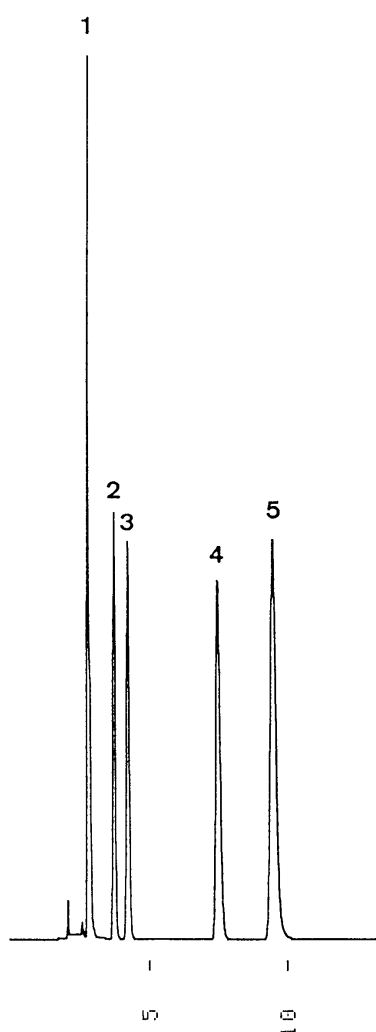


Cefadroxil

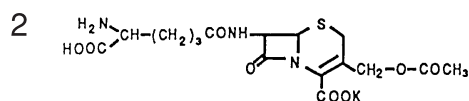


Cephalosporin C potassium

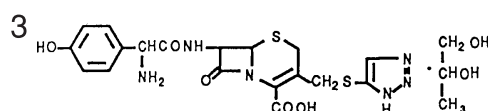

Cefatrizine propylene glycol

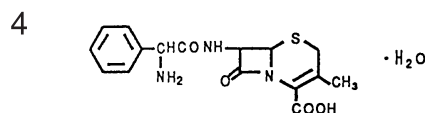

Cephalexin hydrate

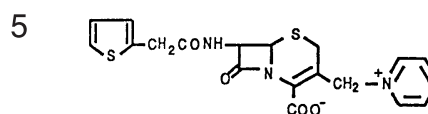

Cephaloridine

Column : J'sphere ODS-L80 (4 μm, 8nm)  
150 × 4.6 mm I.D.  
Eluent : methanol/water/acetic acid (10/85/5)  
Flow rate : 1.0 mL/min  
Temperature : 37 °C  
Detection : UV at 260 nm, 0.16 AUFS  
Injection : 6 μL (0.2 mg/mL)