

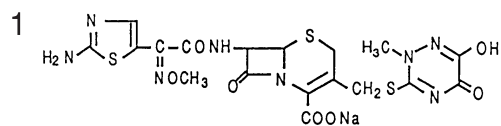


セフェム系抗生物質
Cephalosporin antibiotics

O940929K

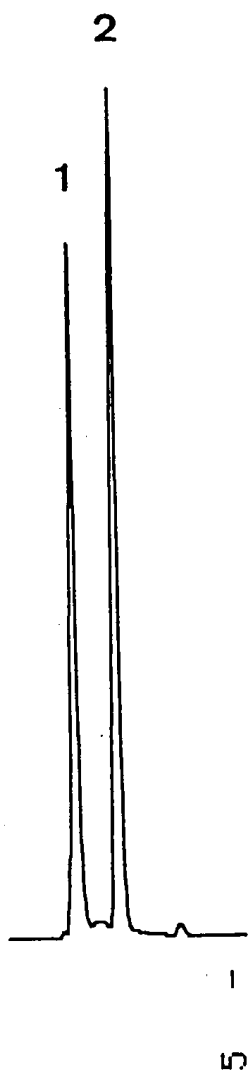

Ceftriaxone sodium

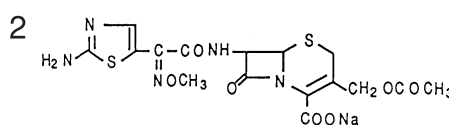

Cefotaxime sodium

Column : J'sphere ODS-H80 (4 μ m, 8nm)
75 \times 2.0 mm I.D.
Eluent : methanol/50mM NH₄H₂PO₄ (30/70)
Flow rate : 0.2 mL/min
Temperature : 37 $^{\circ}$ C
Detection : UV at 260 nm, 0.16 AUFS
Injection : 3 μ L (0.08, 0.16 mg/mL)