


トリアジン系除草剤
Triazine herbicides

O941102F

1. Cyanazine
2. Simazine
3. Atraton
4. Atrazine
5. Prometon
6. Terbumeton
7. Sebuthylazine
8. Propazine
9. Terbutylazine
10. Prometryn
11. Terbutryn

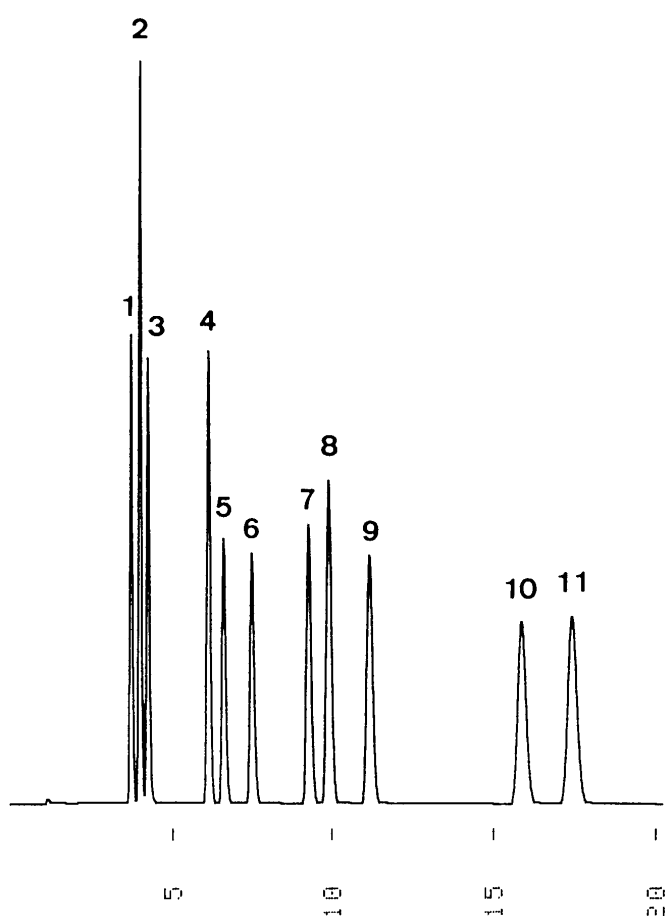

Column : J'sphere ODS-H80 (4 μ m, 8nm)
150 \times 4.6 mm I.D.
Eluent : acetonitrile/50mM NH₄H₂PO₄ (45/55)
Flow rate : 1.0 mL/min
Temperature : 37 $^{\circ}$ C
Detection : UV at 230 nm, 0.32 AUFS
Injection : 5 μ L (0.05 mg/mL)