

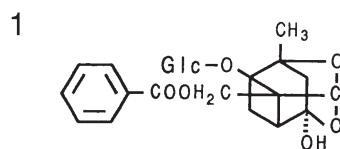


ペオニフロリン
Paeoniflorin

O941121M

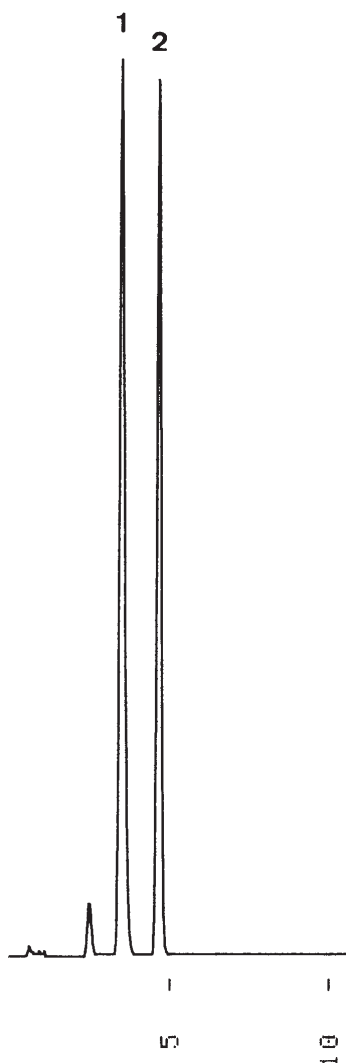

Paeoniflorin

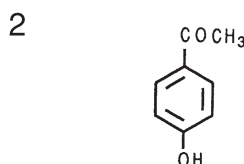

p-Hydroxyacetophenone(I.S.)

Column	: J'sphere ODS-H80 (4 μ m, 8nm) 75 \times 4.6 mm I.D.
Eluent	: THF/50mM (NH ₄) ₂ HPO ₄ (10/90)
Flow rate	: 1.0 mL/min
Temperature	: 37 $^{\circ}$ C
Detection	: UV at 230 nm, 0.16 AUFS
Injection	: 3 μ L (0.67, 0.33 mg/mL)