
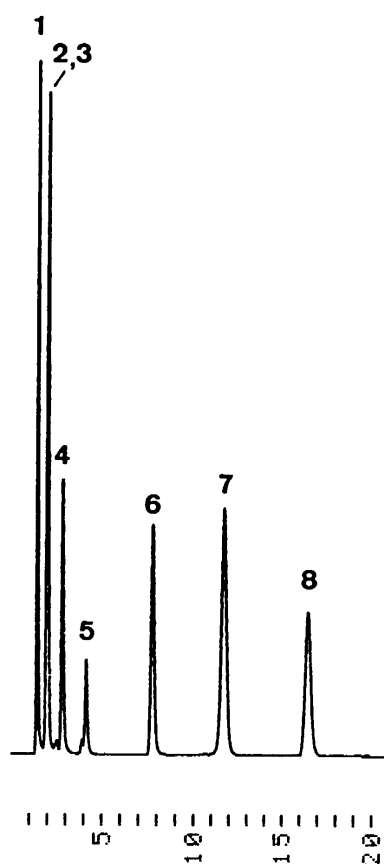
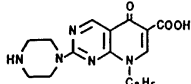
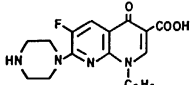
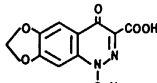
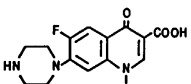
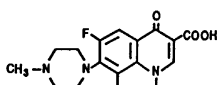
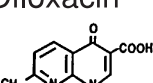
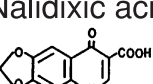


ナリジクス酸系抗生物質
Nalidixic acid antibiotics

O950428A

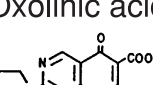

- 1 
Pipemidic acid
- 2 
Enoxacin
- 3 
Cinoxacin
- 4 
Norfloxacin
- 5 
Ofloxacin
- 6 
Nalidixic acid
- 7 
Oxolinic acid
- 8 
Piromidic acid

Column : J'sphere ODS-H80 (4 μ m, 8nm)
150 \times 4.6 mm I.D.

Eluent : 50mM NH₄H₂PO₄ - (NH₄)₂HPO₄ (pH 6.7)/tetrahydrofuran (90/10)
containing 20mM TBA·Br

Flow rate : 1.0 mL/min

Temperature : ambient (25 °C)

Detection : UV at 260 nm, 0.32 AUFS

Injection : 5 μ L (0.03 mg/mL)