

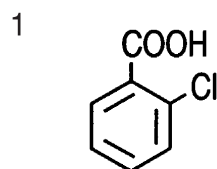


クロロ安息香酸異性体  
Chlorobenzoic acid isomers


O950509B

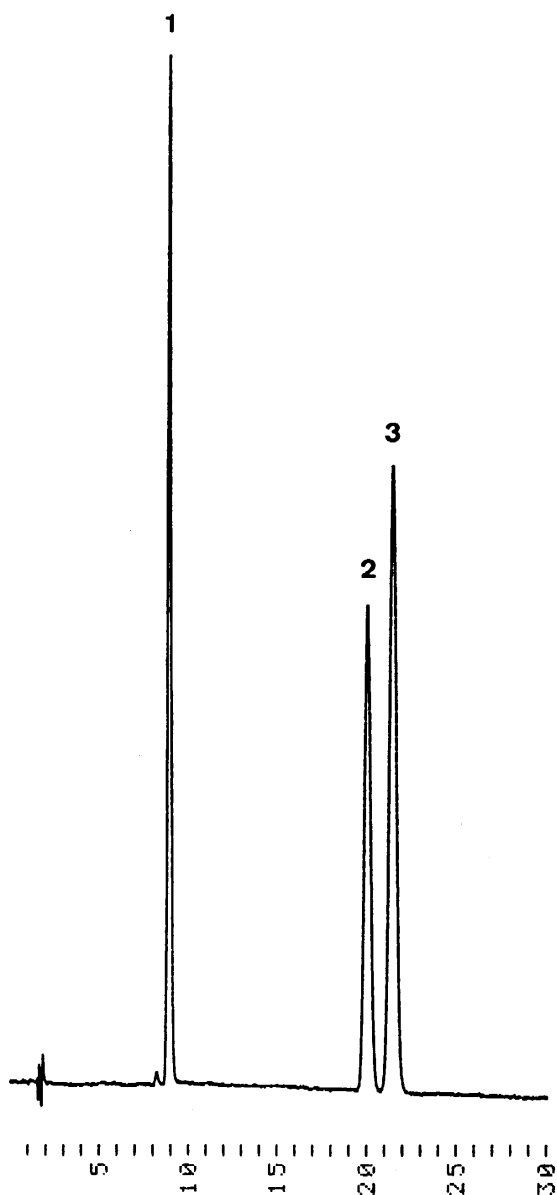

*o*-Chlorobenzoic acid

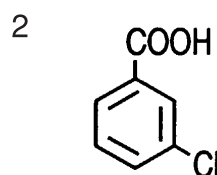

*m*-Chlorobenzoic acid

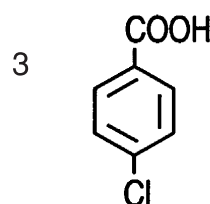

*p*-Chlorobenzoic acid

Column : J'sphere ODS-H80 (4  $\mu$  m, 8nm)  
150  $\times$  4.6 mm I.D.  
Eluent : acetonitrile/water/acetic acid (25/75/1)  
Flow rate : 1.0 mL/min  
Temperature : ambient (27  $^{\circ}$ C)  
Detection : UV at 270 nm, 0.04 AUFS  
Injection : 10  $\mu$ L (0.17 mg/mL)