

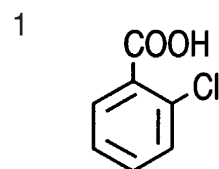


クロロ安息香酸異性体
Chlorobenzoic acid isomers


O950509D

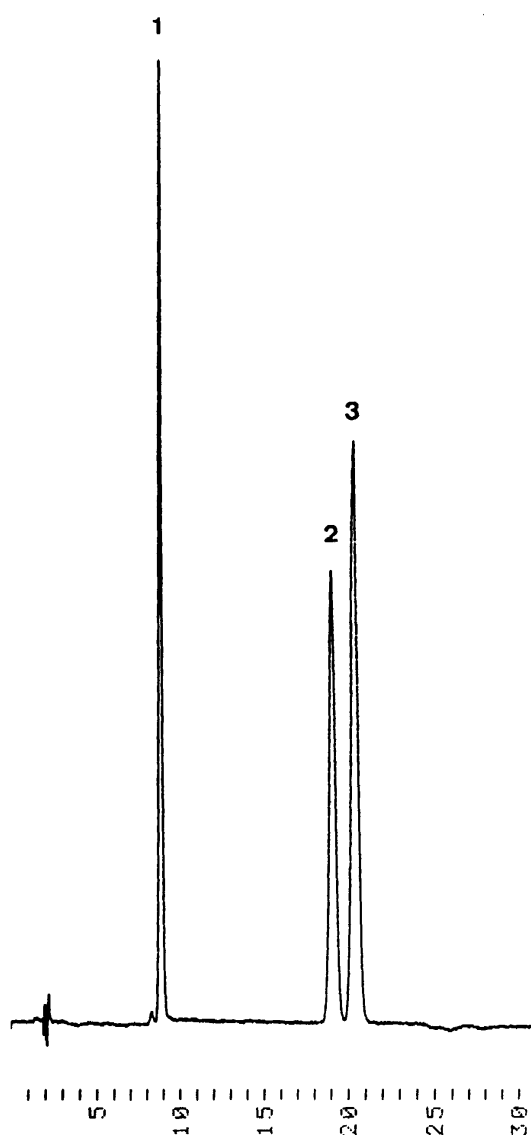

o-Chlorobenzoic acid

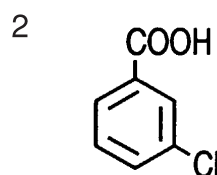

m-Chlorobenzoic acid

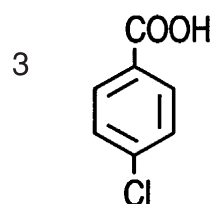

p-Chlorobenzoic acid

Column : YMC-Pack ODS-AM (5 μ m, 12nm)
150 \times 4.6 mm I.D.
Eluent : acetonitrile/water/acetic acid (25/75/1)
Flow rate : 1.0 mL/min
Temperature : ambient (27 $^{\circ}$ C)
Detection : UV at 270 nm, 0.04 AUFS
Injection : 10 μ L (0.17 mg/mL)