
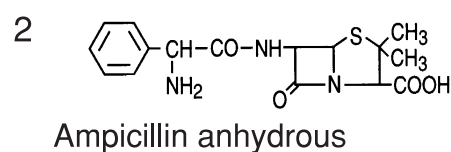
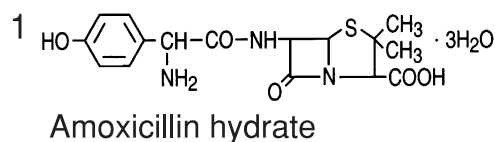


ペニシリン系抗生物質
Penicillin antibiotics

O950920C

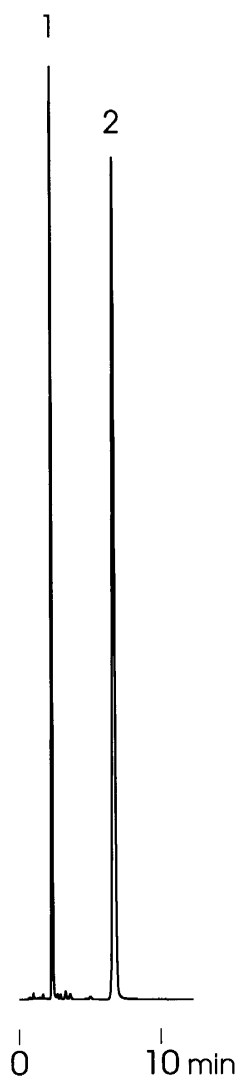

Column : J'sphere ODS-H80 (4 μ m, 8nm)
75 \times 4.6 mm I.D.
Eluent : 50mM NaH₂PO₄ - H₃PO₄ (pH 2.3)/acetonitrile (80/20)
containing 20mM sodium 1-heptanesulfonate
Flow rate : 1.0 mL/min
Temperature : 37 $^{\circ}$ C
Detection : UV at 220 nm, 0.13 AUFS
Injection : 2 μ L (0.3, 0.7 mg/mL)