

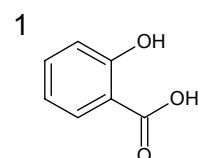


芳香族カルボン酸
Aromatic carboxylic acids


Q090805B

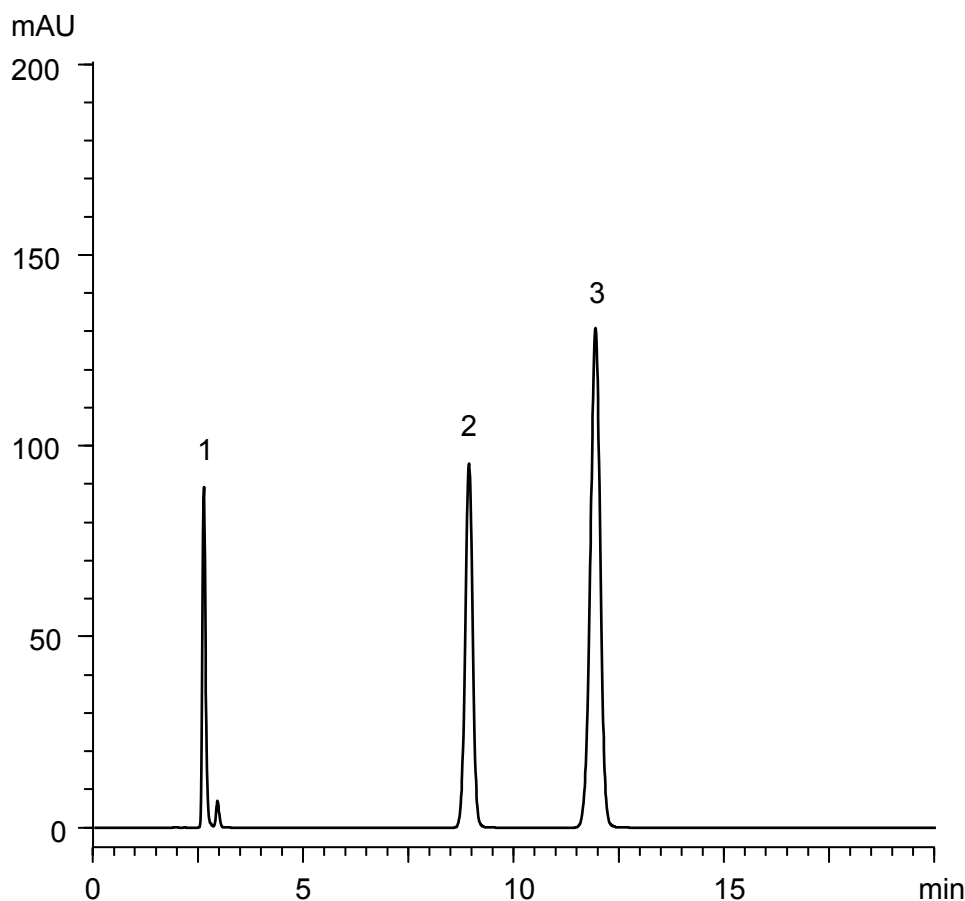

Salicylic acid

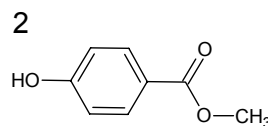

Methyl *p*-hydroxybenzoate (I.S.)

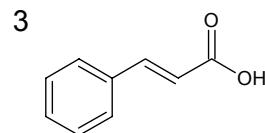

trans-Cinnamic acid

Column	: YMC-Triart C18 (5 μm, 12 nm) 150 X 3.0 mm I.D.
Eluent	: 10 mM CH ₃ COOH-CH ₃ COONH ₄ (pH 4.2)/acetonitrile (75/25)
Flow rate	: 0.425 mL/min
Temperature	: 40°C
Detection	: UV at 254 nm
Injection	: 4 μL (0.02 ~ 0.3 mg/mL)