




ビピリジニウム系除草剤(ジクワット、パラコート)  
Bipyridinium herbicides (diquat, paraquat)

R091201B

1.

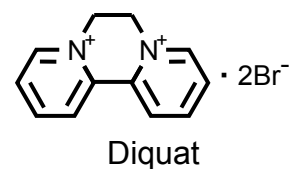

2.

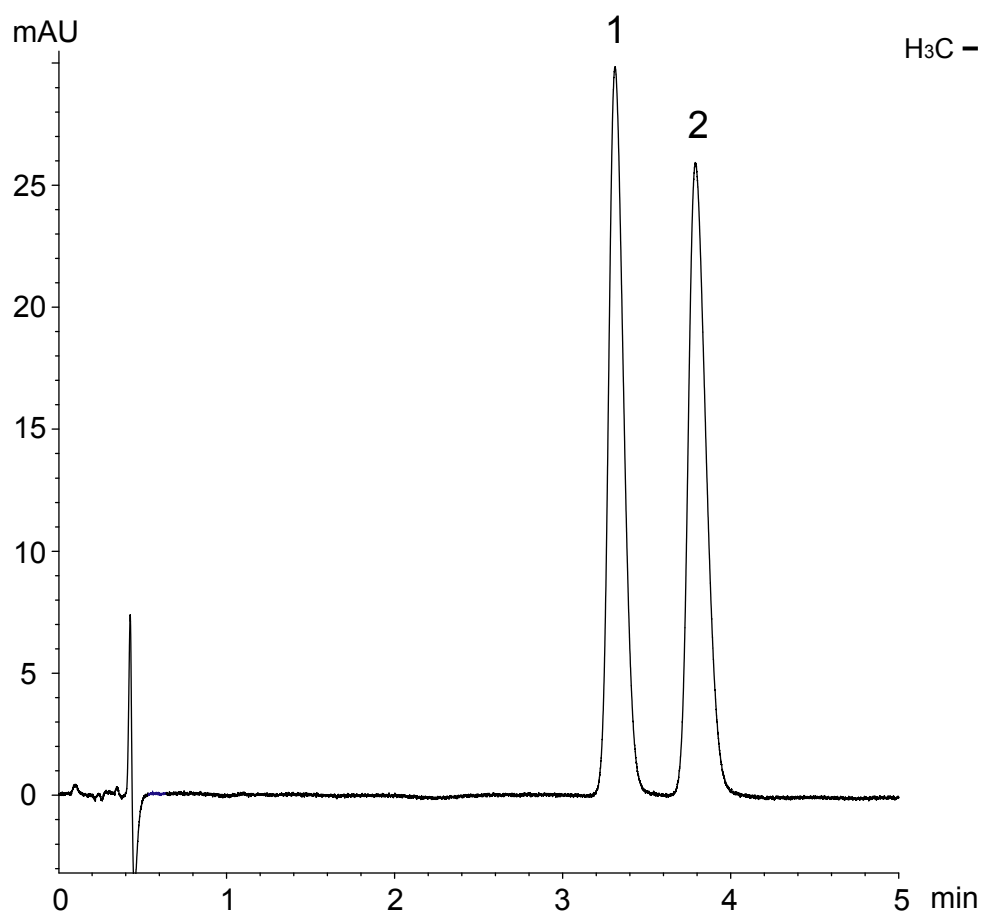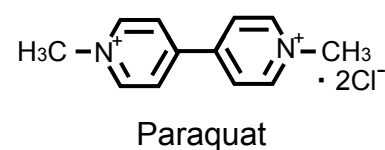

Column : YMC-Triart C18 (3  $\mu$ m, 12 nm)  
50 X 2.0 mmI.D.  
Eluent : water/acetonitrile/HFBA\* (95/5/0.1)  
Flow rate : 0.4 mL/min  
Temperature : 37°C  
Detection : UV at 290 nm  
Injection : 1  $\mu$ L (0.1 mg/mL)

\*heptafluorobutyric acid