

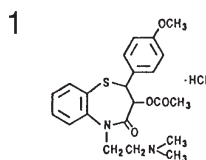


カルシウム拮抗薬  
Calcium antagonists


S930701D

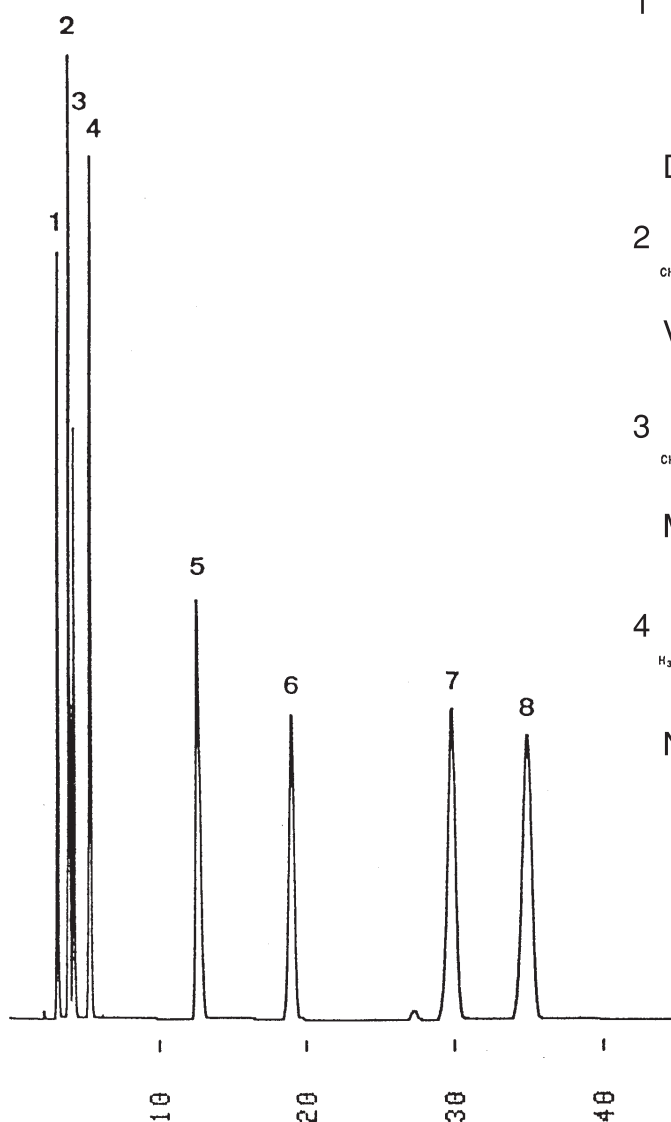


Diltiazem HCl

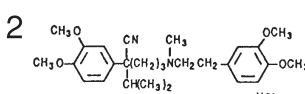


Verapamil HCl

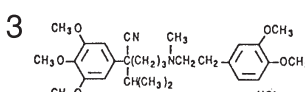


Methoxyverapamil HCl

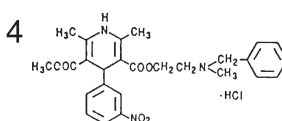


Nicardipine HCl

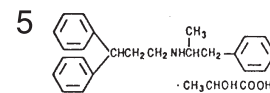


Prenylamine lactate

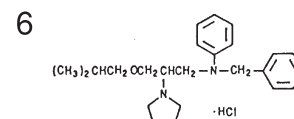

Bepridil HCl

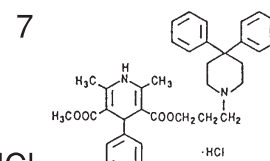

Niguldipine HCl

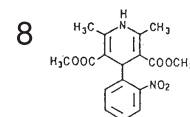

Nifedipine

Column : J'sphere ODS-M80 (4  $\mu$  m, 8nm)  
150  $\times$  4.6 mm I.D.  
Eluent : 20mM KH<sub>2</sub>PO<sub>4</sub> - H<sub>3</sub>PO<sub>4</sub> (pH 3.0)/acetonitrile/THF/TBA·Br  
(70/20/10/0.35, v/v/v/w)  
Flow rate : 0.8 mL/min  
Temperature : 30 °C  
Detection : UV at 270 nm, 0.064 AUFS  
Injection : 8  $\mu$  L (0.05~1.4 mg/mL)