


エリスロマイシン  
Erythromycin antibiotic drugs

S930910B

- 1.Erythromycin
- 2.Erythromycin ethylsuccinate
- 3.Erythromycin estolate

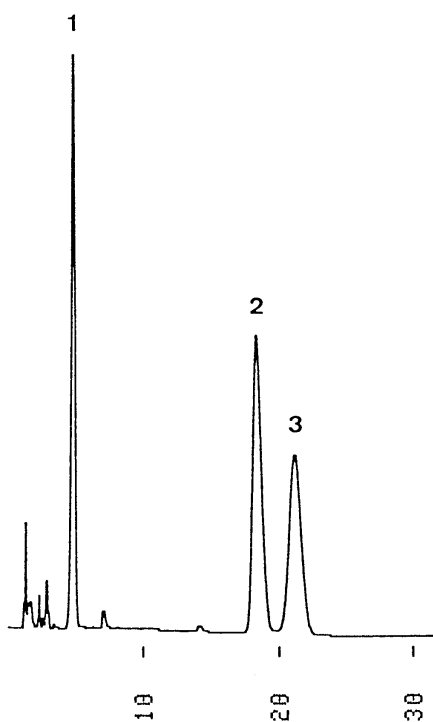

Column : J'sphere ODS-H80 (4  $\mu$ m, 8nm)  
150  $\times$  4.6 mm I.D.  
Eluent : acetonitrile/methanol/100mM (NH<sub>4</sub>)<sub>2</sub>HPO<sub>4</sub> (45/15/40)  
Flow rate : 1.0 mL/min  
Temperature : 45 °C  
Detection : UV at 210 nm, 0.26 AUFS  
Injection : 15  $\mu$ L (2.5 mg/mL)