
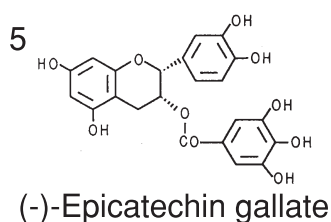
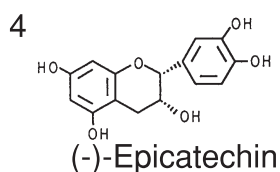
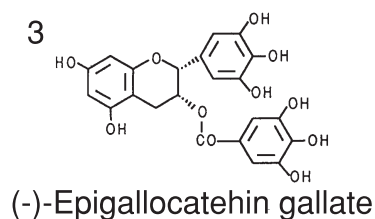
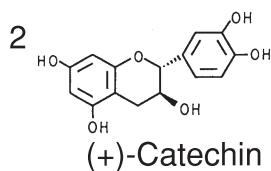
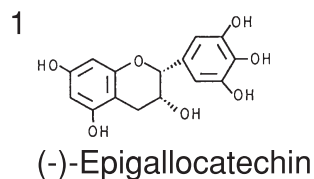


カテキン類  
Catechins

S940124A

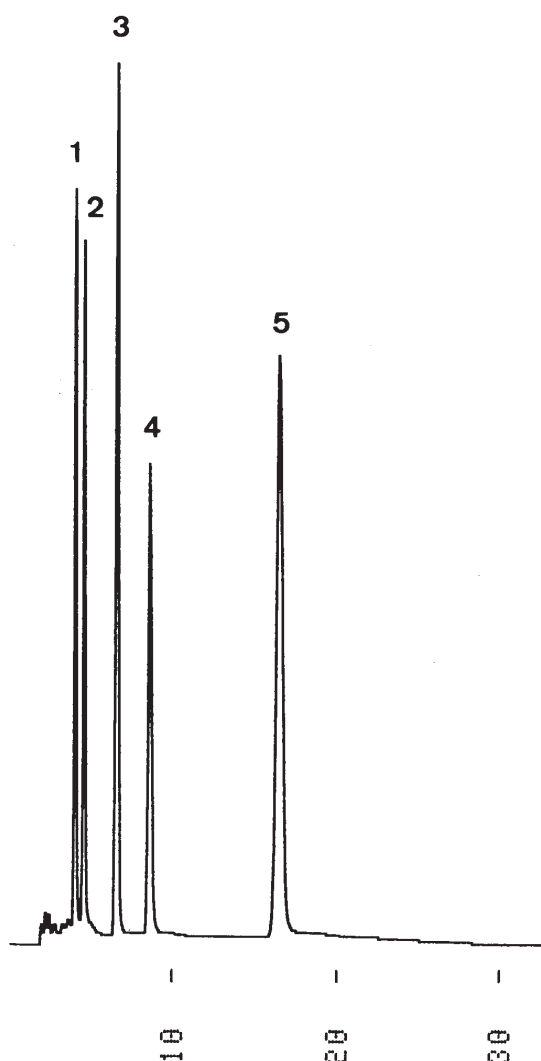

Column : J'sphere ODS-H80 (4  $\mu$ m, 8nm)  
150  $\times$  4.6 mm I.D.  
Eluent : 20mM KH<sub>2</sub>PO<sub>4</sub> - H<sub>3</sub>PO<sub>4</sub> (pH 2.4)/methanol (75/25)  
Flow rate : 0.8 mL/min  
Temperature : 37  $^{\circ}$ C  
Detection : UV at 280 nm, 0.13 AUFS  
Injection : 8  $\mu$ L (0.1~0.4 mg/mL)