


マルトオリゴ糖(G11~G15)
Maltooligosaccharides(G11~G15)

T910718A

1. Maltoundecaose (G11)
2. Maltododecaose (G12)
3. Maltotridecaose (G13)
4. Maltotetradecaose (G14)
5. Maltopentadecaose (G15)

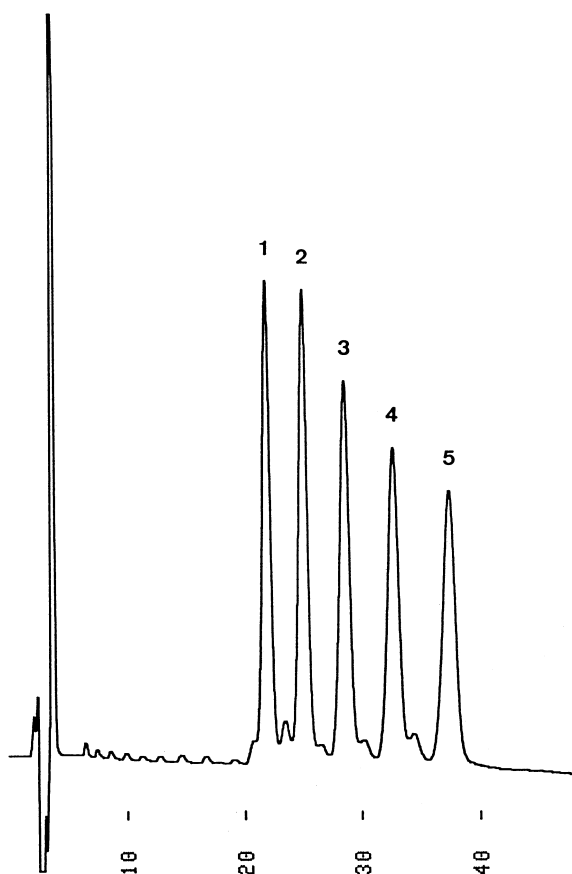

Column : YMC-Pack Polyamine II
250 × 4.6 mm I.D.

Eluent : acetonitrile/water (55/45)

Flow rate : 1.0 mL/min

Temperature : 26 °C

Detection : RI, 32 × 10⁻⁶ RIU/FS

Injection : 200 μL (0.71~0.95 mg/mL)

Sample : Manufactured by Nakano Vinegar Co., Ltd.,
Nakano Central Research Institute