





グリメピリド
Glimepiride

U101117B

A) System suitability solution*
(0.0025 mg/mL *cis*-Glimepiride, 0.5 mg/mL Glimepiride)

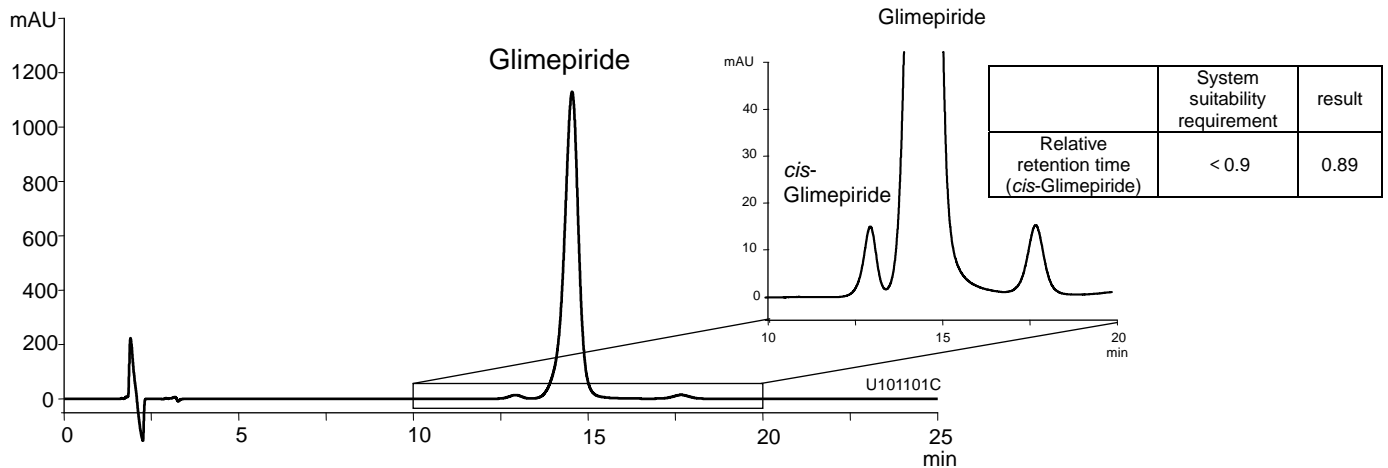

B) Test solution*
(0.5 mg/mL Glimepiride)

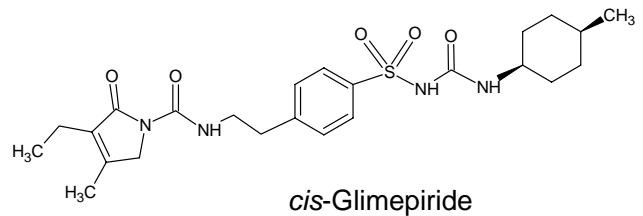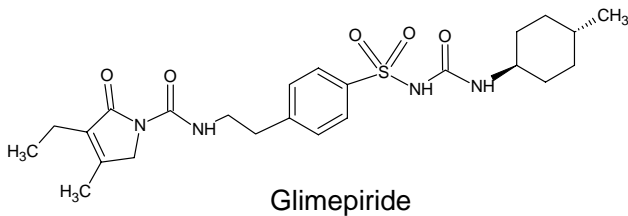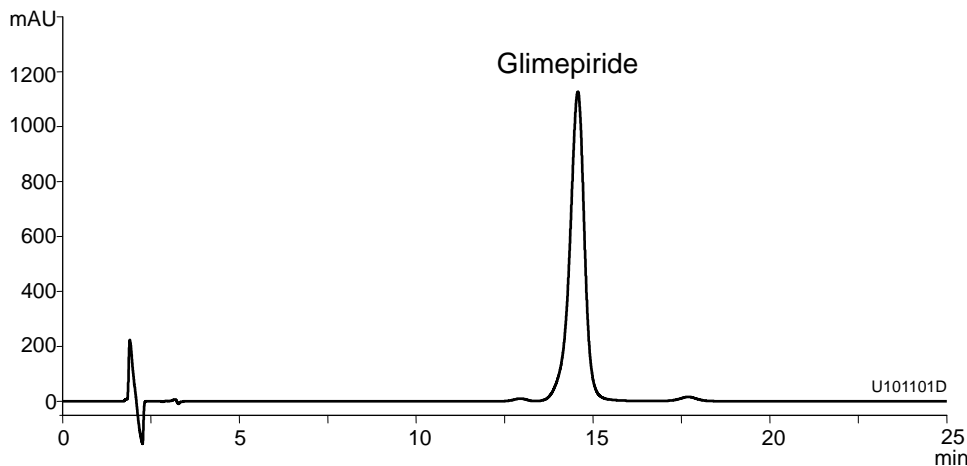

Column : YMC-Pack Diol-NP (5 μm, 12 nm)
150 X 3.0 mmI.D.
Eluent : hexane/2-propanol/acetic acid (900/100/1)
Flow rate : 0.5 mL/min
Temperature : 25
Detection : UV at 228 nm
Injection : 10 μL
(The United States Pharmacopoeia 34th ; *cis*-isomer)

* All system suitability and test solutions were prepared from Glimepiride supplied as a reagent for laboratory use.