




パロキセチン錠  
Paroxetine tablets

U110215A

A) Standard solution\*  
(0.1 mg/mL Paroxetine hydrochloride)

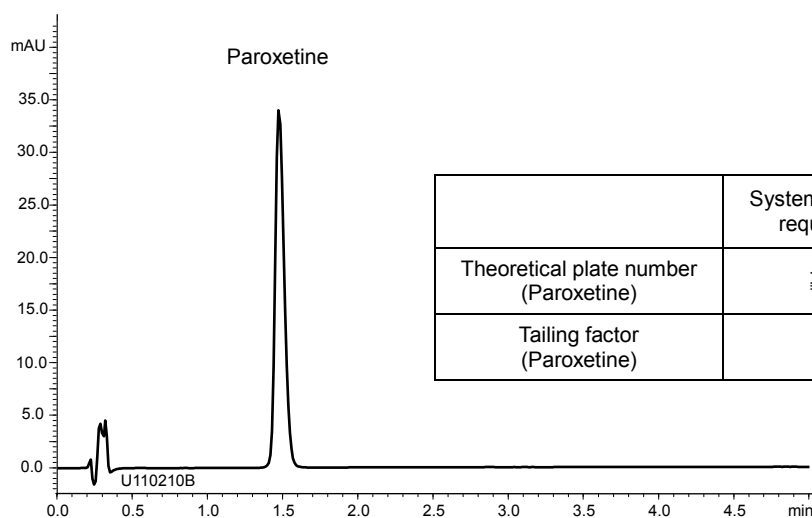

B) Assay solution\*  
(0.1 mg/mL Paroxetine hydrochloride)

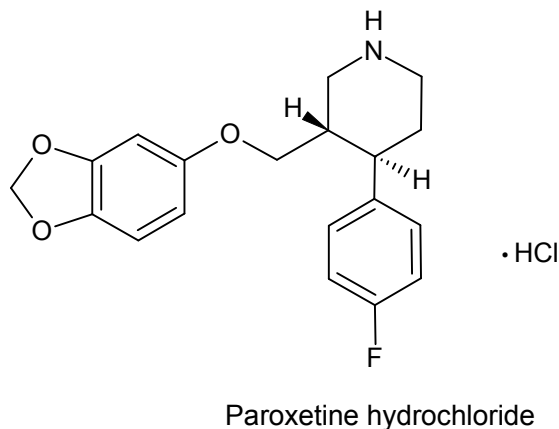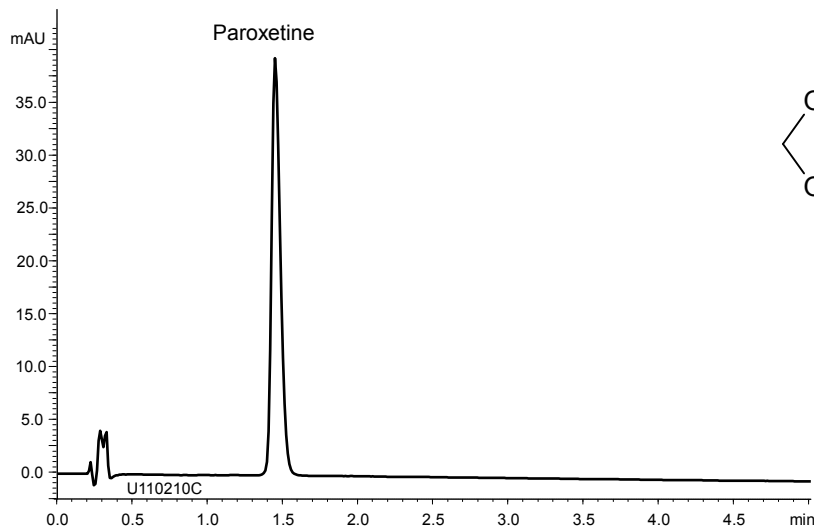

Column : YMC-Pack Pro C8 (3  $\mu$ m, 12 nm)  
35 X 4.6 mm I.D.  
Eluent : acetonitrile/buffer (30/70)  
buffer: water/phosphoric acid/triethylamine (100/0.6/0.3)  
Flow rate : 2.0 mL/min  
Temperature : 25°C  
Detection : UV at 295 nm  
Injection : 5  $\mu$ L  
(The United States Pharmacopoeia 32nd)

\*Standard solution was prepared from Paroxetine HCl supplied as a reagent for laboratory use.  
Assay solution was prepared from Paroxetine HCl tablets.