


モノクローナル抗体  
Monoclonal antibody (MAb)

U170221C3

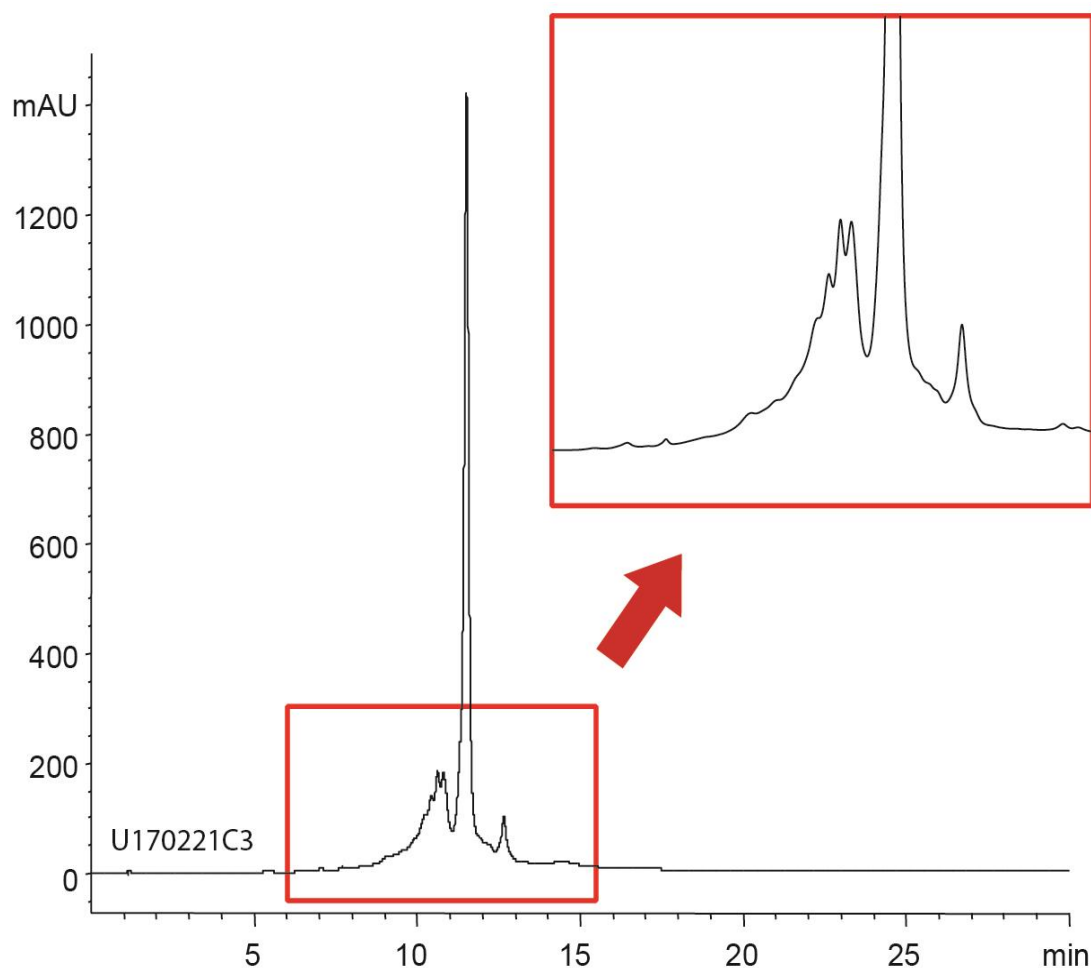

Column	: BioPro IEX SF (3 $\mu$ m) 100 X 4.6 mm I.D.
Eluent	: A) 20 mM NaH <sub>2</sub> PO <sub>4</sub> -Na <sub>2</sub> HPO <sub>4</sub> (pH 6.8) B) 20 mM NaH <sub>2</sub> PO <sub>4</sub> -Na <sub>2</sub> HPO <sub>4</sub> (pH 6.8) containing 0.2 M NaCl 0-50%B (0-30 min)
Flow rate	: 0.5 mL/min (180 cm/hr)
Temperature	: 25°C
Detection	: UV at 215 nm
Injection	: 10 $\mu$ L (2.5 mg/mL Humanized monoclonal IgG1)