



ステロイド  
Steroids

X150731E

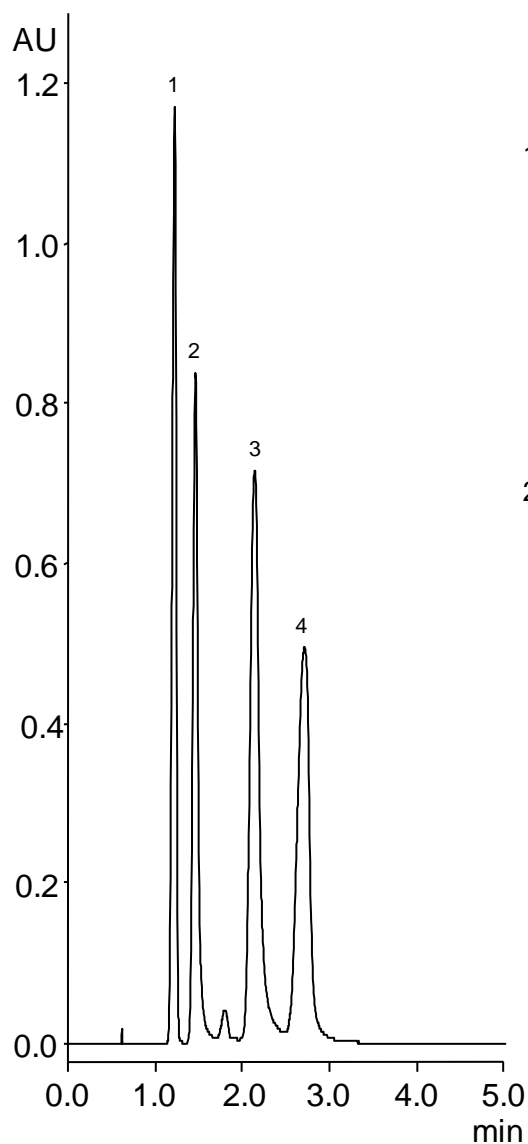


1.

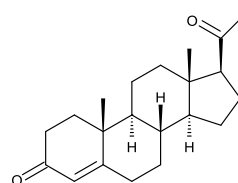


Progesterone

3.

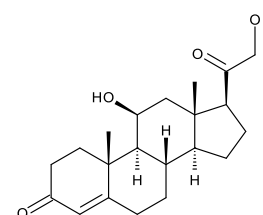


Corticosterone

2.

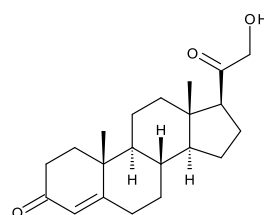

Deoxycorticosterone

4.

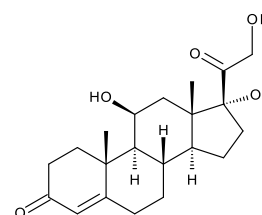

Hydrocortisone

Column	: Alcyon SFC CN (5 μm, 12 nm) 150 X 4.6 mm I.D.
Eluent	: CO <sub>2</sub> /methanol (90/10)
Flow rate	: 3.0 mL/min
Temperature	: 40°C
Detection	: UV at 254 nm
Back Pressure	: 13.8 MPa (2000 psi)
Injection	: 5 μL (0.8 mg/mL)